AWARD NUMBER: W81XWH-13-1-0073

TITLE: 7T Magnetization Transfer and Chemical Exchange Saturation Transfer MRI of Cortical Gray Matter: Can We Detect Neurochemical and Macromolecular Abnormalities?

PRINCIPAL INVESTIGATOR: Seth A. Smith

CONTRACTING ORGANIZATION: Vanderbilt University

REPORT DATE: October 2014

TYPE OF REPORT: Annual

PREPARED FOR: U.S. Army Medical Research and Materiel Command

Fort Detrick, Maryland 21702-5012

DISTRIBUTION STATEMENT: Approved for Public Release;
Distribution Unlimited

The views, opinions and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy or decision unless so designated by other documentation.

REPORT DOCUMENTATION PAGE

Form Approved OMB No. 0704-0188

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Department of Defense, Washington Headquarters Services, Directorate for Information Operations and Reports (0704-0188), 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to any penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number. PLEASE DO NOT RETURN YOUR FORM TO THE ABOVE ADDRESS.

| 1. REPORT DATE | 2. REPORT TYPE | 3. DATES COVERED |
|---------------------------------------|----------------------------------|--|
| 10/21/2014 | Annual | 30 Sep 2013 - 29 Sep 2014 |
| 4. TITLE AND SUBTITLE | | 5a. CONTRACT NUMBER |
| 7T Magnetization Transfer a | and Chemical Exchange Saturation | |
| Transfer MRI of Cortical Gr | ay Matter: Can We Detect Neuro- | 5b. GRANT NUMBER - W81XWH-13-1-0073 |
| chemical and Macromolecular | Abnormalities? | |
| | | 5c. PROGRAM ELEMENT NUMBER |
| | | |
| 6. AUTHOR(S) | | 5d. PROJECT NUMBER |
| Seth Smith | | |
| | | 5e. TASK NUMBER |
| | | |
| | | 5f. WORK UNIT NUMBER |
| E-Mail: seth.smith@vanderbilt.edu | | |
| 7. PERFORMING ORGANIZATION NAME(S | i) AND ADDRESS(ES) | 8. PERFORMING ORGANIZATION REPORT NUMBER |
| The Vanderbilt University | | NOMBER |
| Nashville, TN 37240-0001 | | |
| Nashville, in 3/240-0001 | | |
| | | |
| | | |
| 9. SPONSORING / MONITORING AGENCY | NAME(C) AND ADDDECC/EC) | 40 SPONSOR/MONITOR'S ACRONYM/S) |
| 9. SPONSORING/ MONITORING AGENCY | NAME(3) AND ADDRESS(ES) | 10. SPONSOR/MONITOR'S ACRONYM(S) |
| U.S. Army Medical Research and Ma | ateriel Command | |
| Fort Detrick, Maryland 21702-5012 | atonor communa | 11. SPONSOR/MONITOR'S REPORT |
| Tort Detrick, Maryland 21702-3012 | | NUMBER(S) |
| | | |
| 12 DISTRIBUTION / AVAIL ABILITY STATE | MENT | |

Approved for Public Release; Distribution Unlimited

13. SUPPLEMENTARY NOTES

14. ABSTRACT

We present the first annual report for this project. We have developed and deployed a quantitative MRI set of MRI scans at 7T in healthy volunteers and have developed an analysis pipeline for the quantification of macromolecular and metabolic indices reflective of demyelination and neurotransmitter/protein accumulation. All quantitative MRI methods deployed in humans were optimized and vetted via phantom work and simulation. We are happy to report initial findings that show excellent image quality and differences between the normal cohort and one MS patient. Additionally, we have developed, in conjunction with Dr. Newhouse, a neurocognition expert, an extensive neuropsychiatric battery for correlation studies in year 2. We have had excellent recruitment and retention and look forward to group comparisions in year 2.

15. SUBJECT TERMS -Not Provided

| 16. SECURITY C | LASSIFICATION OF: | | 17. LIMITATION OF ABSTRACT | 18. NUMBER OF PAGES | 19a. NAME OF RESPONSIBLE PERSON USAMRMC |
|----------------|-------------------|--------------|-------------------------------|------------------------|---|
| a. REPORT | b. ABSTRACT | c. THIS PAGE | — | | 19b. TELEPHONE NUMBER (include area code) |
| U | U | U | UU | | |

Table of Contents

| <u>-</u> | <u>Page</u> |
|---|-------------|
| | |
| 1. Introduction | 1 |
| 2. Keywords | 1 |
| 3. Overall Project Summary | 2-9 |
| 4. Key Research Accomplishments | 9-10 |
| 5. Conclusion | 10 |
| 6. Publications, Abstracts, and Presentations | 10 |
| 7. Inventions, Patents and Licenses | 10 |
| 8. Reportable Outcomes | 10-11 |
| 9. Other Achievements | 11 |
| 10. References | 11-12 |
| 11. Appendices | 12+ |

INTRODUCTION:

We recognize that many patients with Multiple Sclerosis (MS) suffer from cognitive impairment at some point in their disease course. However, characterization cognitive change in patients with MS has been difficult to pinpoint, and is hampered by poor quantitative markers. We have two hypotheses: 1) conventional imaging is insensitive to gray matter (GM) changes known to exist in patients with MS, and 2) ultra-high MRI field strengths (7T) would allow an opportunity to study the myelination and metabolic changes of the cortical GM in patients with MS and known cognitive impairment. The purpose of this proposal is to develop and implement a targeted quantitative magnetization transfer (qMT) and chemical exchange saturation transfer (CEST) MRI imaging paradigm at 7T to detect and quantify the level of myelin loss (qMT), protein/peptide changes (amide proton transfer CEST), neurotransmitter deficiencies (GluCEST) in the GM of patients with MS, and to relate these findings to neuropsychiatric evaluation outside the MRI scanner. The scope is to: 1) develop novel, high-resolution, high field, quantitative MRI methods sensitive to myelination and neurochemicals for implementation in the cortical GM of human populations, 2) deploy these methods in patients with MS, 3) relate these findings to measures of cognitive impairment, and 4) develop a lower MRI field strength alternative for direct patient impact.

KEYWORDS:

- Magnetic Resonance Imaging (MRI)
- 7 Tesla (7T)
- Chemical Exchange Saturation Transfer (CEST)
- Magnetization Transfer (MT)
- Brain
- Cortical Gray Matter (cGM)
- Multiple Sclerosis (MS)
- Functional MRI (fMRI)
- Pool Size Ratio (PSR)
- Amide Proton Transfer (APT)
- Glutamate (Glu)
- Myoinositol (mI)
- Cognitive Impairment

OVERALL PROJECT SUMMARY

Task 1. IRB Preparation and Human Subjects Approvals. Completed

Task 2. Develop, optimize and implement advanced, quantitative Magnetization Transfer (MT) and Chemical Exchange Saturation Transfer (CEST) in phantoms and evaluate minimum achievable resolution and the associated reliability of derived indices

Simulation/Phantom Studies

The objective of this task was to develop a best-practice MT and CEST acquisition scheme to be deployed in Task 3 in healthy controls.

Summary of Results/Progress and Accomplishments

CEST -

APT-CEST – We have optimized through simulation and phantom studies a single-power, whole brain APT CEST acquisition for deployment in healthy volunteers and patients (Y2Q1). We began with the protocol presented by Jones et al (1) for whole brain coverage at 7T and increased the in-plane resolution, the coverage,

and modified the saturation scheme, readout, and fat-saturation (now a binomial excitation pulse) to minimize distortions in the phantoms while maintaining sufficient contrast to noise (CNR) for the APT CEST signatures in a reasonable scan time (9 minutes). The relevant scan parameters are as follows:

- Whole Brain (33 slices) 3D Gradient Echo (FFE) with multi-shot EPI (factor = 7) readout
- 1.5 x 1.5 x 5mm³ acquired resolution
- CEST RF Saturation B1 = 2μ T x 25ms (each)
- 64 offset frequencies (Dw = -5 5ppm, $\Delta \omega$ step = 0.2ppm + 14 no Saturation Scans)
- Total Scan Time = 9:10

The simulations were performed using a 3 pool model of the CEST effect presented by Zaiss et al (2) and inputting estimates for T1, T2, pool sizes and exchange rates for the macromolecular, labile (amide protons) protons, and water as follows: Simulations were designed to identify optimal CEST preparation (RF irradiation power and bandwidth) parameters for APT contrast. All simulations were carried out utilizing the scripting environment in MATLAB 2012b (Mathworks, Natick, Massachusetts) on an Apple iMac (Cupertino, CA; 3.0 GHz, dual core CPU). Theoretical saturation was modeled according to the Bloch equations for three pools: bulk water (free), semisolid macromolecular (conventional MT), and mobile macromolecular (CEST) pools. This was achieved using the simple matrix solution to numerically solve the Bloch equations (3). This model assumes a T2-dependent Super-Lorentzian absorption lineshape for the macromolecules (4). Physical values of exchange rate, T1, and T2, were fixed according to (5.6). The bulk water was modeled as a Lorentizan with T1/T2 = 1538 ms/45ms. The semisolid macromolecular pool was modeled with T1/T2/exchange rate = 1600ms/0.01ms/20 Hz with offset = -2.34ppm (7). The labile proton pool was modeled as a Lorentzian with T1/T2/exchange rate = 500ms/20ms/50Hz with offset = 3.5ppm and a concentration of 0.001 compared to bulk water (1.0) and macromolecular pool (0.1). Therefore our simulations included a single RF irradiation, a brief delay for spoiling, an on-resonance excitation, and a delay for readout. The B₁ amplitude (power) was varied over amplitudes from 1 μT to 3 μT while holding the duration constant at 25 ms. The B_1 amplitude was subsequently fixed to 1 µT while the pulse duration was varied from 0 ms to 60 ms.

For phantom scans, we performed the above paradigm at various in-plane resolutions as low as 0.75mm², but determined through curve analysis that at these resolutions the CNR was insufficient to parse out the CEST effect from the background noise. From these simulations and phantom studies, the above pulse sequence paradigm was chosen to maximize CEST contrast derived from APT.

GluCEST – For glutamate, we performed the same simulations as above, but we modeled the off-resonance saturation as given in (8), and assumed an exchange rate for glutamate = 100 Hz, pool size = 0.001 compared to bulk water (1.0) and macromolecular concentration (0.01) and $\Delta\omega$ = 3.0ppm. For the GluCEST acquisition, due the concern of overlapping resonances (GABA, Glutamate, and other Amines) a high-spectral resolution acquisition needed to be obtained ($\Delta\omega$ = 0.2ppm spacing), thus scan time becomes prohibitive for extremely high resolution. However, it should be pointed out that the amine resonances that may reside in juxtacortical WM will be significantly less than the adjacent GM, so a slightly poorer resolution acquisition will not be problematic if a high spectral resolution scan is obtained. We therefore, decided to utilize a scan very similar to that which as been presented by Dr. Reddy (8,9) and thus we will implement a single slice GluCEST acquisition in vivo at a resolution of 1.9 x 1.9 x 5mm³.

- Single Slice (2D) Gradient Echo (FFE) with multi-shot TFE (40 shots) readout
- 1.9 x 1.9 x 5mm³ acquired resolution
- CEST RF Saturation B1 = 4.25μ T x 10ms (each) x 100segments at 90% duty cycle
- 50 offset frequencies ($\Delta \omega = -5 5$ ppm, $\Delta \omega$ step = 0.2ppm)
- Total Scan Time = 11:36

qMT -

We have chosen to implement the selective inversion recovery (SIR) quantitative MT (qMT) to quantitatively extract the pool size ratio (PSR), which has been shown to be reflective of myelin. We have developed the SIR

approach at 7T as discussed in the Q1Y1 and Q2Y1 progress reports. This pulse sequence has been shown in previous reports, published (10) and provided in Appendix 3. However, in Q3Y1, we studied via simulation the impact of partial volume effects where we know there is a non-negligible MT effect in GM and certainly a strong effect in WM. Thus, in juxtacortical voxels where a blend of GM and WM may occur, poorer in-plane resolution results in an inability for the model to remain stable when deriving the PSR values (i.e. two different PSR values may fit equally well when there are two populations within a voxel). Therefore, we proposed a reduced number of slices but increased the in-plane resolution. From simulations, we feel that this provides the most robust acquisition method to be deployed in patients. Thus, from our phantom studies, we have determined that a 1 x 1 x 2mm³ acquisition with 5 slices sampled at 14 TIs (TI = 6ms, 10ms, 16ms, 26ms, 42ms, 68ms, 110ms, 178ms, 288ms, 468ms, 760ms, 1233ms, 2000ms, 8000ms) will be performed in patients with MS and healthy controls.

- Inversion prepared 3D Gradient Echo (FFE) with multi-shot TFE (2 shots) readout
- 1 x 1 x 2mm3 acquired resolution (5 slices)
- 14 Inversion Times (TI = 6ms, 10ms, 16ms, 26ms, 42ms, 68ms, 110ms, 178ms, 288ms, 468ms, 760ms, 1233ms, 2000ms, 8000ms) at a constant delay time (TD = 2500ms)
- Total Scan Time = 10:11

Conclusion of Task 2: Simulation and Phantom-optimized qMT and CEST acquisitions

Through simulation and phantom studies, we have devised a final protocol to be deployed in healthy volunteers and patients with MS. A summary of the protocol is given below, and a complete protocol is given in Appendix 1.

Final Summary of Protocol implemented in healthy controls (Task 3) and patients with MS (Task 5).

- Constant RF APT CEST 9:10
- Constant RF GluCEST 11:36
- SIR qMT 10:11
- Bloch-Siegert B1 mapping 1:42
- Dual-echo B0 mapping :04
- T1w MPRAGE Anatomical 2:12
- fMRI Resting State 8:34
- fMRI N-Back task 8:30
- fMRI Trailmaking task 4:14

The current scan time for all scans is approximately 1 hour.

Task 3 – Implement current best practice for MT and CEST in healthy volunteers and evaluate reliability The objective of this task was to implement a best-practice MT and CEST acquisition scheme in healthy controls.

Summary of Results/Progress and Accomplishments

The above protocol has been implemented in 20 healthy volunteers at the close of year 1. We have additionally repeated this paradigm in 8 healthy volunteers. We have further scheduled the remainder of the healthy controls to be scanned in the coming month. There is one delay to report in that at the close of Year 1, our SOW stated that we would have recruited 50 healthy volunteers into the study. As pointed out in quarterly reports Q1Y1 and Q2Y1, we struggled initially making the phantoms to study the impact of resolution on the final protocol. This resulted in less than the expected 50 healthy controls. However, we have already scheduled these remainder healthy volunteers, and have a 40% return rate on for repeat visits to understand the reproducibility. We will complete the healthy volunteers and repeat visits in Q1Y1

In a follow-up to the Q3Y1 report, we have added three fMRI scans in collaboration with Dr. Paul Newhouse, our neurocognition expert. It is important that we note that this does not change the scope, but rather offers a unique opportunity to study the cognitive function in the MRI in healthy volunteers and patients for greater understanding of the relationship between the advanced, quantitative measures and outside-scanner neurocognitive battery. This is exceptionally unique as neither of these three fMRI scans have been studied in MS patients with known cognitive impairment and creates an exceptionally rich data set to mine for understanding neurocognitive decline in MS patients.

In 20 healthy volunteers (with 8 repeat acquisitions), we have obtained the entire proposed MRI protocol as given above. Preliminary results follow under Task 4.


Task 4 – Analyze the derived indices in healthy volunteers and evaluate reproducibility (1 month) The objective of this task was to 1) develop an analysis pipeline for constructing maps and deriving indices reflective of GM and juxtacortical WM from the quantitative MRI acquisitions prepared in Task 3, and 2) to ascertain these indices in preparation for analysis of reproducibility.

Summary of Results/Progress and Accomplishments

Develop an analysis pipeline for routine analysis of data generated.

As the first part of Task 4, we well understood the need to 1) correct patient motion in an individual scan (motion-correction), 2) co-register data across scans into the same space for robust analysis (co-registration), and 3) segmentation of WM and GM for histogram and descriptive statistics of each derived index.

WM/GM segmentation was performed in FAST using 3 classes as implemented in the FSL toolbox (FMRIB, Oxford, UK). The co-registration was performed using FNIRT (non-linear registration, FSL, FMRIB, Oxford,

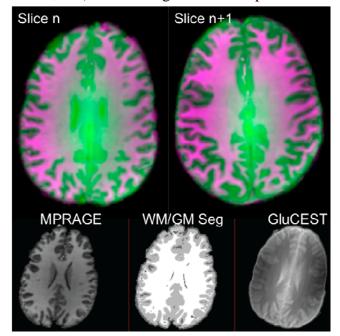

UK) to put the quantitative maps into the MPRAGE space such that the segmentation can be applied. To that end, Figure 1, shows two slices of an APT-CEST acquisition motion corrected and co-registered. The magenta color indicates the WM, and the green indicates the GM derived from the MPRAGE anatomical acquisition and overlaid on the APT map. It can clearly be seen that the agreement between the MPRAGE and the APT maps is high and the WM and GM clearly seen. The bottom panels show the process of joining the MPRAGE, the WM/GM Segmentation and the GluCEST-weighted acquisition for completeness.

Figure 1 – (Top) coregistration and segmentation results applied to APT maps, and (bottom) the target (MPRAGE), segmented map, and GluCEST-weighted acquisition.

APT CEST analysis and results

We have constructed the APT CEST maps in the following manner. First, the CEST spectrum for each voxel is normalized, corrected for B1 drift and fit to a single-lorentizian (11) and the minimum spectral intensity is shifted to an offset ($\Delta\omega$) = 0 for B0 correction. After this correction, the difference between the data and the fit create a Lorentzian residual. The residual between $\Delta\omega$ = 3.25 and 3.75 are integrated and termed the APT Lorentzian. An example of this is shown in Figure 2A left panel. To assess reproducibility we created a single

CEST spectrum for the whole brain and compared visit 1 to visit 2. Figure 2A right panel shows the results of two healthy repeated studies, where the green/blue curves are subject 1, and magenta/black curves are subject two over two times points. As it can be seen, the reproducibility is high over the whole brain. Individual structure assessments are ongoing. It can also be appreciated the spectral quality of the CEST spectrum using this analysis approach.

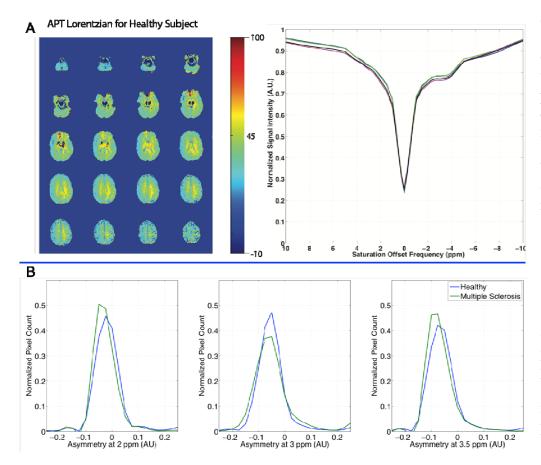



Figure 2: (A) APT maps for all slices derived in a healthy volunteer and concomitant test-retest CEST spectra in 2 healthy volunteers. The test-retest is over the whole brain, and the green/blue and magenta/black spectra pairs are from the same volunteers. (B) histogram analysis of segmented GM in 20 healthy volunteers and 1 MS patient with clinically noted cognitive impairment. Note that histograms at 2ppm and 3.5ppm show a downward shift of the MS patient relative to the healthy control indicating initial sensitivity to the pathology of cortical GM damage.


Once the data were segmented, we compared the GM averaged over all healthy volunteers (n = 20) and 1 MS patient with clinically diagnosed cognitive impairment histogram via analysis shown in Figure 2B. We performed the same calculation, but also examined resonances at 2ppm (hydroxyl and sensitive to myoinositol) and 3ppm (amine protons sensitive to **GABA** and Glutamate) and 3.5ppm transfer (amide proton sensitive to pH and protein concentration). It should be pointed out that the 1 MS patient examined here is actually part of Task 5, but it is important to show here as it points towards the sensitivity of the measurement. In this patient, there is obvious downward shift of the GM histograms at 2ppm and 3.5ppm giving the impression that we are detecting cortical and perhaps even some subcortical changes in protein concentration and myoinositol. In Figure 2B, the green is the MS patient, and the blue is the

average over healthy volunteers. We are exceptionally excited by this initial result and felt it important to share here as Task 5 will indeed prove the sensitivity of these advanced techniques to MS.

Glutamate CEST (GluCEST) analysis and results

GluCEST analysis proceeded as presented in (8,9). We performed GluCEST analysis in 20 healthy volunteers and 4 MS patients at the time of this report, though only one had been analyzed and is presented here. In short, GluCEST-weighted images were collected for a single slice with high spectral saturation fidelity (see Task 2) at a slice slightly superior to the corpus callosum. Sample GluCEST weighted images are shown as a function of offset frequency in Figure 3, top panel. From these maps, the GluCEST spectra were corrected for B0 and B1 in the fashion presented in (9), and a GluCEST asymmetry map at $\Delta\omega = 3.0$ ppm was generated. Figure 3, bottom panel shows two healthy volunteers and one patient with MS, clinically diagnosed with cognitive

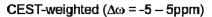

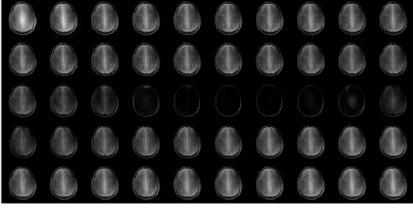

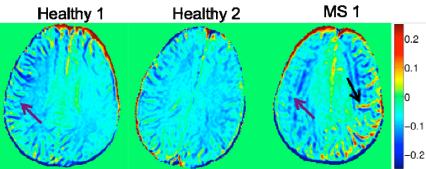


Figure 3: (top) example GluCEST-weighted images as a function of offset frequency (bottom) comparison of GluCEST maps for two healthy volunteers and one MS patient. Note the apparent differences between the MS patient and healthy volunteers (black and magenta arrows)

impairment. First, it can be seen that the GluCEST maps show excellent contrast between WM and GM, with the GM having higher GluCEST signal than the WM (expected). What is exciting to note (and in general, part of Task 5) is the visual differences between the patient and the two controls shown here. The MS patient shows elevated GluCEST signal (black arrow) on the left side, but apparently diminished GluCEST signal in the right cortical GM (magenta arrow) compared to healthy subject 1 (magenta arrow). This seems to indicate, at least at the early stages, that GluCEST is detecting cortical GM differences between healthy and MS patients.

As with the APT CEST, we examined the entire 20 healthy control cohort in comparison to the MS case and Figure 4A shows the average GluCEST spectra derived from GM and WM in healthy patients (blue and black, respectively) compared to the segmented GM in the MS case (red). It can be seen that the spectral quality is high and there is visual difference between the spectra for healthy and MS GM. Further, we analyzed the histogram of GluCEST signals for all GM

voxels in all healthy volunteers and compared that to the 1 MS patient clinically diagnosed with cognitive challenges (Figure 4B). As with the APT CEST, it can be seen that the MS patient shows a downward trend compared to the healthy volunteer indicating the possibility of being sensitive to cortical GM pathology, which

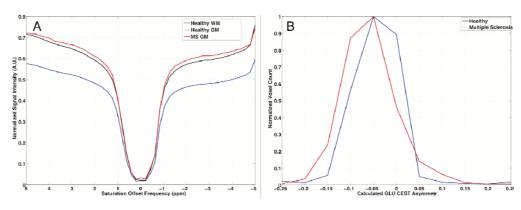


Figure 4: (A) Average GluCEST spectra for WM and GM in healthy (blue black) and MS patient (red) GM. (B) Histogram analysis shows that the MS has a downward shift of the GluCEST signal compared to the healthy volunteer (blue)

we will study in detail in Task 5. When we started this project, we decided against an "all-in-one" CEST acquisition scheme and rather have deployed two CEST acquisitions. One sensitive to APT (and apparently myoinositol) and one sensitive to glutamate. From Figure 4B, and in comparison to Figure 2B, the histograms from the APT-CEST analysis at $\Delta \omega = 3.0$ ppm (glutamate) show now difference between healthy controls and the 1 MS patient, however, when

utilizing the GluCEST acquisition, in the same patient, a difference can be appreciated. The rationale for this is

that the sensitivity to exchanging species is determined by the power of the RF CEST saturation. For the APT, the exchange rate is on the order of 20-100Hz, whereas for glutamate amines, the exchange rate is faster (50-200Hz). Thus, to be maximally sensitive to both, two separate pulse powers are necessary. We discovered this as part of Task 2 in the phantoms and are proud to note that it was the right choice going forward.

Quantitative Magnetization Transfer (qMT) analysis and results

High-resolution selective inversion recovery (SIR) qMT was performed in the same cohort as for CEST and analyzed according to (10) and given in Appendix 3. In short, an inversion recovery MRI sequence was performed using a modified inversion pulse that is relatively insensitive to B1 and B0 inhomogeneities. The inversion times were selected to sample the bi-exponential recovery known to exist when magnetization transfer is present. For every voxel, the SIR signal equation was fit to the recovery curve and the exchange rate (kmf), pool size ratio (PSR) and longitudinal relaxation time (R1f) was fit. Appendix 3 provides the manuscript that contains details of the pulse sequence, and fitting method.

Here we report the initial analyses and results from the newly deployed method. Figure 5, left panel shows the PSR, R1f, and kmf maps from a representative healthy volunteer and one patient with MS and concomitant cognitive impairment. From the PSR maps, it can be seen the high level of discrimination between WM (yellow/orange) and GM (blue) which agrees well with the R1f maps. Two things should be noted for the MS patient. First, the WM shows a globally decreased PSR which is indicative of demyelination across the entire slice, while the R1f and kmf maps do not show a similar pattern (discussed next). Secondly, when looking at the GM, it is not apparent that the patient and the control have different PSR values, yet MS patients are known

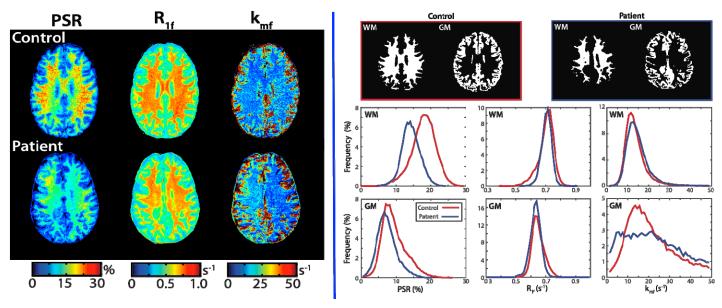


Figure 5: (left) qMT-derived maps for a healthy control (top) and MS patient (bottom). (right) segmentation results (top) and histogram analyses for the average healthy controls (red) and the MS patient (blue). Note the downward shift of the PSR for both WM and GM in this patient, while R1 and kmf are indistinguishable.


to have myelin loss in the GM. However, Figure 5, right panel, shows a histogram analysis of over the healthy volunteers (N = 20) and 1 MS patient. The top row shows the segmented WM and GM for an example healthy volunteer, and the bottom panels show the histograms for WM and GM for all of the qMT-derived indices. It can be seen again that for WM, the MS patient has a substantially downshifted PSR, normal R1f and kmf.

Importantly, however, in the GM, the patient also shows a small downshift of the PSR, while R1f and kmf are not markedly different (it should be noted that kmf showed some instabilities in this patient). This is important

to note in that one argument about MT imaging is that it is hypothesized that R1 drives the change in the MT effect moreso than does the macromolecular content. This figure shows that rather, in both WM and GM, the PSR is abnormal, but neither R1 nor kmf are indicating the sensitivity to WM and GM macromolecular pathology. Task 5 will explore this further when we examine a larger MS cohort.

fMRI analysis and results

In addition to the quantitative measures that have already been shown and at the advice of our mentor, Dr. Paul Newhouse, a neurocognition expert, we added 3 fMRI acquisitions to the MRI paradigm. Those three methods were an N-Back task, a Trail-making task, and a resting-state fMRI acquisition. The N-back and Trail-making tasks are important as they are also performed outside the MRI scanner, so we will be able to provide direct correlations between what is performed in the MRI and outside of the MRI. This further allows us to directly,

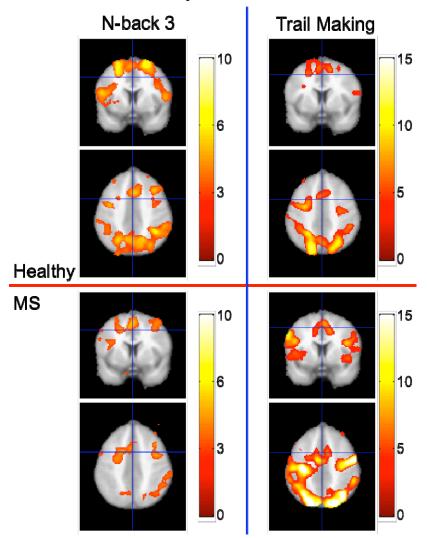

and non-invasively probe cognitive performance in a manner that is not only unique, but it has not been performed in the MS population at 7T. We are encouraged by the initial results and wish to present those here. We performed the fMRI in 20 healthy volunteers and 1 MS patient at the time of this report. Figure 6 shows a direct comparison of the N-back 3 (left panels) and the Trail-making (right panels) for a single healthy volunteer (top panels) and a patient with MS (bottom panels). As it can be seen, at the same significance threshold, there are activation differences between the healthy volunteer and the MS patient. These are especially noted for areas in the superior cortex where working memory is targeted. A greater confidence will be gained with a group analysis, but this will be reserved for Task 5.

Figure 6 (left) N-back 3 and (right) Trail-making fMRI activation patterns in a healthy (top panels) volunteer and one patient with MS (bottom panels). Note for the N-back, there is less activation in the MS patient, where as for the Trail-making task, there are greater areas of activation.

Discussion of Task 4

We are pleased with the quality of data that has been generated and are actively enrolling healthy volunteer sand patients with MS. We have scanned 20 healthy volunteers, and 4 MS patients to this point, and have had 8 healthy volunteers return for a 2nd visit. We are slightly behind the enrollment expectation, but this is in large part due to the ground work necessary to start the human studies. We are, however, encouraged that because of this extensive focus on sequence optimization, that the data quality remains superb in all subjects. We have scheduled the remainder of the healthy subjects, developed a robust pipeline for analysis, and have shown initial success in implementing these acquisitions in patients with MS. We expect no delays in finalizing enrollment and will continue with reproducibility analysis from which we can gauge the level of expected deviation from normal in patients with MS.

Task 5 – Implementation in Patients with MS and concomitant cognitive impairment

The objective of this task is to deploy, and analyze the MRI acquisitions in patients with MS.

We have implemented the MRI paradigm in 4 patients with MS and expect no delays in enrollment for the remainder of the MS patients. The details of preliminary results are given in Task 4 with the goal of comparing both image quality, and quantitative differences. We have nothing else to report for Task 5 at the time of this annual report.

Task 6 – Cross-sectional analysis of derived indices between patients with MS and healthy volunteers and correlation with clinical measures of cognitive impairment derived from the Minimal Assessment of Cognitive Function in MS (MACFIMS)

The objective of this task is to compare quantitative MRI indices across cohorts, implement neuropsychiatric evaluations in both healthy and MS patient cohorts, and derive correlations with quantitative MR indices.

Neuropsychiatric Assessment Battery (Outside MRI Scanner)

In collaboration with Dr. Paul Newhouse, we have decided to additionally obtain neuropsychiatric evaluations in healthy volunteers in addition to patients with MS. This will provide us with a greater understanding of the variance across cohorts. Therefore, we have obtained neuropsychiatric data using the paradigm below in 20 healthy volunteers (8 repeats) and 4 MS patients.

Tasks - Outside the scanner and BEFORE coming to Vanderbilt

Questionnaire and survey already developed in REDCap to be completed at home and in a calm environment. These surveys will collect data related to baseline mood, anxiety, and cognitive profile.

Tasks - Outside the scanner at Vanderbilt (< 1 hour total)

- Short measure of day-of mood/anxiety
- N-back test (2-back or 3-back): measures working memory
- PASAT: measures working memory
- Trail making test (both A and B): measures planning/executive function
- "Black Box" (choice reaction time, critical flicker fusion; pre-scan and post-scan): measures processing speed/reaction time
- *Buschke selective reminding test (8 trials): measures include encoding and long-term memory
- Digit Symbol Substitution Test/DDST: measures visual memory
- Posner cueing task: measures attentional shift

KEY RESEARCH ACCOMPLISHMENTS:

- 1. Developed a Chemical Exchange Saturation Transfer (CEST) simulation pipeline to model the effects of glutamate, amide proton transfer, and myoinositol in the gray matter (GM) at 7T. These simulations reflect contributions from the metabolite of interest, the magnetization transfer (MT) effect, and the direct water saturation.
- 2. Developed a quantitative magnetization transfer (qMT) simulation pipeline to model the effects of macromolecular concentration of the white matter (WM) and GM in healthy tissue and tissue impacted by multiple sclerosis (MS). This pipeline incorporates modeling of the semisolid fraction while taking into account relaxation times (T1 and T2) changes that are known to occur in each cohort.
- 3. Phantoms have been created that attempt to model the in vivo scenario. That is, they have varied concentration of metabolite and concomitant concentration of semi-solid components (MT) at constant pH. These phantoms include: glutamate, myoinositol, glycogen, bovine serum albumin (BSA), and agarose. These phantoms can be leveraged for greater understanding of the in vivo results

- 4. Developed and optimized a set of novel MRI acquisition strategies to study the CEST effects of glutamate, myoinositol, amide protons, taking into account corrections for B1 and B0 inhomogeneities. These have been deployed in vivo.
- 5. Developed and optimized a high-resolution (1mm² in-plane resolution) qMT acquisition that is of sufficient resolution to assess cortical GM and juxtacortical WM in healthy and MS cohorts.
- 6. Developed a multi-modal motion-correction, coregistration and WM/GM/CSF segmentation strategy that not only maps the confidence of the measurements made in small structures, but also puts all acquisitions in the same space for descriptive statistics on each cohort.
- 7. Implemented a set of fMRI experiments to assess working memory and resting-state fluctuations in patients as compared to healthy controls. This will work in conjunction with item #8: neuropsychiatric evaluation
- 8. Created a detailed neuropsychiatric evaluation paradigm for assessing cognitive impairment that can be performed outside of the scanner environment and can be leveraged for correlation testing in the final year of this project.
- 9. In summary, we have developed, optimized, and deployed a multi-parametric, multi-modal MRI toolkit to assess neurochemical, macromolecular, functional, and structural changes in vivo at 7T. Additionally, we have developed a detailed neuropsychiatric evaluation paradigm to utilize for comparisons across cohorts, which can be further extended to any MRI study of cognitive impairment.

CONCLUSION:

We have concluded the first year of this project and have many significant contributions to report. First, we have for the first time, developed a battery of quantitative MRI methods that are of sufficient resolution and sensitivity to characterize cortical gray matter in healthy volunteers and patients with multiple sclerosis. To this end, we have a < 1 hour exam card that can be deployed on any 7T scanner that can investigate neurochemical composition, macromolecular/myelin deficiency, and functional impairment in a patient cohort. Additionally, we have generated data that suggests that there are differences between healthy volunteers and patients with MS, while expected, has yet to be shown because lower MRI field strengths have insufficient sensitivities to these macromolecules and neurochemicals, and insufficient resolution to study only gray matter. We have additionally partnered with a neurocognitive expert, and, with his help, developed a novel neuropsychiatric battery to assess cognition in MS. We understand that these techniques, while not novel, have not been implemented in patients with MS, and may provide evidence for greater scope in any patient with neurocognitive decline. We have studied 20 healthy volunteers and a handful of MS patients to this point and will expand the patient and control enrollment in year 2. We have finally, developed a pipeline for analysis that requires minimal human interaction and will deploy this for real-time analysis in year 2. It is important to note, however, that while these techniques are currently being explored for use at 7T, in year 2, we will develop a lower field, and thus, significantly more clinically relevant, set of exams that will provide a similar MRI toolkit. Lastly, the exam as developed here is not specific for MS and can be implemented in a wide range of patients and volunteers to explore the neurochemistry, functional processing, and macromolecular composition of cortical gray matter.

PUBLICATIONS, ABSTRACTS, AND PRESENTATIONS:

- 1. Lav Press: Nothing to report
- 2. Peer-Reviewed Scientific Journals: Nothing to report
- 3. **Invited Articles**: Nothing to report
- 4. **Abstracts**: 7 abstracts to the International Society for Magnetic Resonance in Medicine (ISMRM) annual conference have been prepared and will be submitted in Q1Y2 (November 12, 2014 deadline)

INVENTIONS, PATENTS AND LICENSES:

REPORTABLE OUTCOMES:

- 1. High-resolution, optimized 7T MRI acquisition strategies (so-called Exam Card in Philips language) to quantitatively evaluate the macromolecular, metabolic, functional, and structural characteristics that can (and currently is) be implemented in healthy controls and any patient cohort. A complete listing of the MRI acquisition paradigm (Exam Card) is given in Appendix 1.
- 2. Assessment of the reproducibility and stability of each measurement over time is currently ongoing.
- 3. A complete set of neuropsychiatric assessments, some of which are completely novel in patients with MS
- 4. Analysis pipeline for generation of quantitative MRI-derived indices. The pipeline includes motion correction, multi-modal image co-registration, and WM/GM/CSF segmentation along with generation of quantitative indices for further statistical comparisons.
- 5. A CEST simulation GUI for further studies of the CEST effect in vivo.
- 6. Collection of experiments, simulations, and phantom studies that have provided evidence for the minimal resolution attainable while maximizing sensitivity to change in patient populations.
- 7. Complete set of MRI and neuropsychiatric data in 20 healthy volunteers (with 8 additional repeat visits) and 4 patients with MS and clinically confirmed cognitive impairment. These data are summarized in the Overall Project Summary.

OTHER ACHIEVEMENTS:

- 1. Because of the nature of the experiments performed as a part of this grant, that is to implement the highest resolution quantitative MRI at 7T in WM and GM of healthy participants and MS patients, we have been able to extend these tools to the spinal cord, which has significant impact and scope for patients with other forms of neurological injury. One manuscript (*Dula AN, Pawate S, Dethrage LM, Conrad BN, Barry RL, Smith SA. CEST of the Cervical Spinal Cord at 7 Tesla. Submitted to NMR in Biomed on 30-Sept-2014*) has already been submitted on the results from this extension to other parts of the body.
- 2. The phantoms and simulations that were created are not MRI field strength dependent. Therefore, we have changed the simulations to study the impact of transitioning to a lower-field strength (i.e. 3T) for greater clinical implementation. We have begun to study the sensitivity for 3T utilization. This is critical, as it was noted in our initial application, that 7T MRI scanners are not directly clinically impactful. Thus, we have been able to create an MRI acquisition strategy at lower field strength which will be deployed in year 2. We are excited about the possibility of reaching a greater clinical community with the studies in year 2 at 3T.
- 3. From the preliminary results generated in the first year of this award, Dr. Pawate (co-investigator) is preparing a grant submission to the National Multiple Sclerosis Society in February 2015 to study cerebral changes in primary progressive MS (PPMS) patients.

REFERENCES: List all references pertinent to the report using a standard journal format (i.e. format used in Science, Military Medicine, etc.)

- 1. Jones CK, Polders D, Hua J, Zhu H, Hoogduin HJ, Zhou J, Luijten P, van Zijl PC. In vivo three-dimensional whole-brain pulsed steady-state chemical exchange saturation transfer at 7 T. Magn Reson Med 2012;67(6):1579-1589.3291747
- 2. Zaiss M, Schmitt B, Bachert P. Quantitative separation of CEST effect from magnetization transfer and spillover effects by Lorentzian-line-fit analysis of z-spectra. J Magn Reson 2011;211(2):149-155
- 3. Woessner DE, Zhang S, Merritt ME, Sherry AD. Numerical solution of the Bloch equations provides insights into the optimum design of PARACEST agents for MRI. Magn Reson Med 2005;53(4):790-799
- 4. Morrison C, Stanisz G, Henkelman RM. Modeling magnetization transfer for biological-like systems using a semi-solid pool with a super-Lorentzian lineshape and dipolar reservoir. Journal of magnetic resonance Series B 1995;108(2):103-113
- 5. Michaeli S, Garwood M, Zhu XH, DelaBarre L, Andersen P, Adriany G, Merkle H, Ugurbil K, Chen W. Proton T2 relaxation study of water, N-acetylaspartate, and creatine in human brain using Hahn and Carr-Purcell spin echoes at 4T and 7T. Magn Reson Med 2002;47(4):629-633
- 6. Rooney WD, Johnson G, Li X, Cohen ER, Kim SG, Ugurbil K, Springer CS, Jr. Magnetic field and tissue dependencies of human brain longitudinal 1H2O relaxation in vivo. Magn Reson Med 2007;57(2):308-318
- 7. Hua J, Jones CK, Blakeley J, Smith SA, van Zijl PC, Zhou J. Quantitative description of the asymmetry in magnetization transfer effects around the water resonance in the human brain. Magn Reson Med 2007;58(4):786-793.3707117
- 8. Cai K, Haris M, Singh A, Kogan F, Greenberg JH, Hariharan H, Detre JA, Reddy R. Magnetic resonance imaging of glutamate. Nature medicine 2012;18(2):302-306.3274604
- 9. Cai K, Singh A, Roalf DR, Nanga RP, Haris M, Hariharan H, Gur R, Reddy R. Mapping glutamate in subcortical brain structures using high-resolution GluCEST MRI. NMR in biomedicine 2013;26(10):1278-1284.3999922
- 10. Dortch RD, Moore J, Li K, Jankiewicz M, Gochberg DF, Hirtle JA, Gore JC, Smith SA. Quantitative magnetization transfer imaging of human brain at 7 T. NeuroImage 2013;64:640-649.3625658
- 11. Jones CK, Huang A, Xu J, Edden RA, Schar M, Hua J, Oskolkov N, Zaca D, Zhou J, McMahon MT, Pillai JJ, van Zijl PC. Nuclear Overhauser enhancement (NOE) imaging in the human brain at 7T. NeuroImage 2013;77:114-124.3848060

APPENDICES

```
Philips MRI Protocol Dump
Created on
10/28/2014 10:40:51 AM
```

```
Comment
   Created by ExamCard_to_XML with inputs: "E:\Export\20141021 CEST fMRI.ExamCard" on system (Vanderbilt University :: 192.168.71.10)
    Software Stream
   3.2.1.0
3.2.1.0

Expand All | Collapse All
(2) (1) (1) (13)

Hospital (2)

B20141021 CEST fMRI (13) 52:37.4

SCOUT SHC32 00:28.7

WIP MTX SENSE 32ch 01:28.9

T1_3D_TFE_iso1.25mm_s2.5s SENSE_Sagittal 02:12.2

CEST_Interspersed_3uT 09:09.6

CEST_Reddy_GluCEST 11:37.0

B1_Reddy_multiAngle 01:42.0

B0_Reddy_multiAngle 00:03.9

GMT High Res 02:38.7

FMRI_RESTINGSTATE 08:34.0

FMRI_TRAILMAKING 04:14.0

FMRI_TRAILMAKING 04:14.0

FMRI_TRAILMAKING 04:14.0

T1_3D_TFE_quantGeo 00:55.7

T2star_multiEcho 01:02.7

(B1860723-4F8F-476e-8075-D42C65706693) (0)
      [B1860723-4F8F-476e-8075-D42C65706693] (0)
```

| Hospital (2) 201 | 41021 CEST fMRI | (13) 52:37.4 SCOU | T SHC32 00:28.7 | | |
|--|---|-------------------------|-----------------|----------------------|--------------|
| INFO PA | GE | GEOMETE | RY | CONTRA | ST |
| Total scan duration | 00:28.7 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | M2D |
| Act. TR/TE (ms) | 7.8 / 4.9 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| ACQ matrix M x P | 256 x 128 | Xmit Coil selection | MTX-Volume- | + ZOOM | no |
| ACQ voxel MPS (mm) | 0.98 / 1.95 / | | T/R | Contrast enhancement | T1 |
| | 10.0 | User def elem sel | no | Acquisition mode | cartesian |
| REC voxel MPS (mm) | 0.98 / 0.98 / | element selection | All | Fast Imaging mode | TFE |
| C (0/) | 10.0 | connection | conn-A | shot mode | multishot |
| Scan percentage (%) | 50 | Coil selection 2 | RX-Intf-2 | TFE factor | 64 |
| TFE shots | 2 | element selection | All | startup echoes | default |
| TFE dur. shot / acq (ms) | | Dual coil | yes | +TFE followup echoes | 0 |
| TFE shot interval (ms) | 1063.369 | Multi coil | no | shot interval | shortest |
| Min. TI delay | 287.8819 | CLEAR | no | profile order | linear |
| Act. WFS (pix) / BW (Hz) | 3.513 / 288.4 | FOV FH (mm) | 250 | Echoes | 1 |
| Min. WFS (pix) / Max. | 1.297 / 781.3 | AP (mm) | 250 | partial echo | no |
| BW (Hz) | 1.277 / 701.3 | stack RL (mm) | 50 | shifted echo | no |
| Min. TR/TE (ms) | 7.8 / 2.8 | Voxel size FH (mm) | 0.9765625 | TE | in-phase |
| RF avg power computed | | AP (mm) | 1.953125 | (ms) | 4.93426 |
| (W) | 11000 | Slice thickness (mm) | 10 | Flip angle (deg) | 15 |
| SAR / head | < 100 % | Recon voxel size (mm) | 0.9765625 | TR | shortest |
| Whole body / level | < 0.1 W/kg / | Fold-over suppression | no | Halfscan | no |
| | normal | Reconstruction matrix | 256 | Water-fat shift | user defined |
| B1 rms | 1.32 uT | SENSE | no | (pixels) | 3.5 |
| PNS / level // VUIIS : | 35 % / normal | k-t BLAST | no | Shim | default |
| dortch : | | Stacks | 3 | mDIXON | no |
| Sound Pressure Level | 26.85212 | current | Α | Fat suppression | no |
| (dB) | | type | parallel | Water suppression | no |
| МОТІО | | slices | 3 | TFE prepulse | invert |
| Cardiac synchronization | no | slice gap | user defined | slice selection | no |
| Heart rate > 250 bpm | no | gap (mm) | 10 | shared | no |
| Respiratory | no | slice orientation | sagittal | delay | user defined |
| compensation | | fold-over direction | AP | (ms) | 800 |
| Navigator respiratory | no | fat shift direction | F | PSIR | no |
| Flow compensation | no | Slice scan order | default | +inv pulse type | +default |
| fMRI echo stabilisation | no | Stack scan order | ascend | MTC | no |
| Motion smoothing | no | Move table per stack | no | T2prep | no |
| NSA | 1 | Stack alignment | no | Research prepulse | no |
| DYN/AN | | Stack display order | no | Diffusion mode | no |
| Angio / Contrast enh. | no | PlanAlign | no | Elastography mode | no |
| Quantitative flow | no | REST slabs | 0 | SAR mode | high |
| Manual start | | Catheter tracking | no | B1 mode | default |
| | no | Interactive positioning | no | SAR Patient data | auto |
| +Abuse dynamic loop | no | Allow table movement | no | PNS mode | low |
| Dynamic study | no | OFFC/A | NG | Gradient mode | default |
| Arterial Spin labeling | no | Stacks | 3 | SofTone mode | no |
| POST/PR Preparation phases | | current | A | | - |
| | auto | Stack Offc. AP | 0 | [| |
| Interactive F0 | no | (P=+mm) | | | |
| SENSE ref. scan | no | RL (L=+mm) | 0 | | |
| SmartPlan survey | no | FH (H=+mm) | 0 | | |
| | no | Ana AD (dea) | 0 | | |
| B0 field map | | Ang. AP (deg) | - | | |
| B1 field map | no | RL (deg) | 0 | | |
| B1 field map MIP/MPR | no no | | | | |
| B1 field map MIP/MPR Images | no no M, no, no, no | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image | no no M, no, no, no | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image Calculated images | no no M, no, no, no M no, no, no, no | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | no no M, no, no, no M no, no, no, no Grey matter | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast | no no M, no, no, no M no, no, no, no Grey matter soft | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode | no no M, no, no, no M no, no, no, no Grey matter soft real time | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode Save raw data | no no M, no, no, no M no, no, no, no Grey matter soft real time no | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode Save raw data Hardcopy protocol | no no M, no, no, no M no, no, no, no Grey matter soft real time no no | RL (deg) | 0 | | |
| B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode Save raw data | no no M, no, no, no M no, no, no, no Grey matter soft real time no | RL (deg) | 0 | | |

| ☐ Hospital (2) ☐ 201 | 41021 CEST fMRI | (13) 52:37.4 WIP N | ITX SENSE 32ch | 01:28.9 | |
|------------------------------|-----------------|--|----------------|----------------------|-----------|
| INFO PA | GE | GEOMETR | RY | CONTRA | ST |
| Total scan duration | 01:28.9 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | 3D |
| Act. TR/TE (ms) | 8.0 / 0.75 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| ACQ matrix M x P | 96 x 75 | Xmit Coil selection | MTX-Volume- | loop order | zy_order |
| ACQ voxel MPS (mm) | 5.52 / 7.07 / | | T/R | + ZOOM | no |
| | 6.00 | User def elem sel | no | Contrast enhancement | T1 |
| REC voxel MPS (mm) | 5.52 / 5.52 / | element selection | All | Acquisition mode | cartesian |
| Scan percentage (%) | 3.00 78.125 | connection | conn-A | Fast Imaging mode | none |
| | 1 | Coil selection 2 | RX-Intf-2 | 3D non-selective | no |
| Packages Act. WFS (pix) / BW | 0.489 / 2071.3 | element selection | All | Echoes | 1 |
| (Hz) | 0.469 / 20/1.3 | Dual coil | yes | partial echo | no |
| Min. WFS (pix) / Max. | 0.486 / 2083.3 | Multi coil | yes | shifted echo | no |
| BW (Hz) | 0.1007 2000.0 | CLEAR | no | TE | shortest |
| RF avg power computed | 0.05085949 | FOV FH (mm) | 530 | Flip angle (deg) | 1 |
| (W) . | | RL (mm) | 530 | TR | shortest |
| SAR / head | < 2 % | stack AP (mm) | 300 | Halfscan | no |
| Whole body / level | 0.0 W/kg / | Voxel size FH (mm) | 5.520833 | Water-fat shift | minimum |
| | normal | RL (mm) | 7.066667 | Shim | default |
| B1 rms | 0.18 uT | AP (mm) | 3 | mDIXON | no |
| PNS / level // VUIIS : | 18 % / normal | Recon voxel size (mm) | 5.520833 | Fat suppression | no |
| dortch : | 27 02005 | Image shutter | yes | Water suppression | no |
| Sound Pressure Level (dB) | 27.92985 | Fold-over suppression | no | MTC | no |
| MOTIO | NI . | Slice oversampling | default | Research prepulse | no |
| Cardiac synchronization | | RF select. FOS | no | Diffusion mode | no |
| Heart rate > 250 bpm | no no | Reconstruction matrix | 96 | Elastography mode | no |
| Respiratory | no | SENSE | no | SAR mode | high |
| compensation | 110 | k-t BLAST | no | B1 mode | default |
| Navigator respiratory | no | Overcontiguous slices | yes | SAR Patient data | auto |
| comp | | Stacks | 2 | PNS mode | low |
| Flow compensation | no | current | A slices | Gradient mode | default |
| fMRI echo stabilisation | no | | 100 | SofTone mode | no |
| NSA | 3 | slice orientation | coronal | | |
| SMART | yes | fold-over direction | RL | | |
| DYN/AN | G | fat shift direction | F 1 | | |
| Angio / Contrast enh. | no | Chunks | | · | |
| Quantitative flow | no | Stacks as packages | no | | |
| Manual start | no | Move table per stack | no | | |
| +Abuse dynamic loop | no | Stack alignment Stack display order | no no | · | |
| Dynamic study | no | PlanAlign | no | | |
| Arterial Spin labeling | no | REST slabs | 0 | } | |
| POST/PR | | Catheter tracking | no | } | |
| Preparation phases | full | Interactive positioning | no | } | |
| Interactive F0 | no | Allow table movement | no | } | |
| SENSE ref. scan | yes | OFFC/AI | | | |
| SmartPlan survey | no | Stacks | 2 | } | |
| B0 field map | no | current | A | | |
| B1 field map | no | Stack Offc. AP | 11.5916 | | |
| MIP/MPR | no | (P=+mm) | | | |
| Images | no, no, no, no | RL (L=+mm) | -3.434925 | [| |
| Autoview image | no | FH (H=+mm) | 4.208135 | | |
| Calculated images | no, no, no, no | Ang. AP (deg) | 0 | | |
| Reference tissue | White matter | RL (deg) | 0 | | |
| Preset window contrast | soft | FH (deg) | 0 | | |
| Reconstruction mode | immediate | | | | |
| Save raw data | no | | | | |
| Hardcopy protocol | no | | | | |
| Ringing filtering | default | | | | |
| Elliptical k-space shutter | derauit | | | | |

| Hospital (2) 2014 | A1021 CEST fMDI | (12) 52:27 4 - T1 2F | TEE ico1 25mm | s 2 Ec SENSE Sagittal Of | 0-10-0 |
|---|---|---------------------------------------|---------------|--------------------------|---------------|
| INFO PAG | | GEOMETR | | CONTRA | |
| Total scan duration | 02:12.2 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | 3D |
| Act. TR/TE (ms) | 2.8 / 1.32 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| ACQ matrix M x P | 204 x 204 | Xmit Coil selection | MTX-Volume- | + ZOOM | no |
| ACQ voxel MPS (mm) | 1.25 / 1.25 / | ATTIL COIL SELECTION | T/R | Contrast enhancement | T1 |
| ACQ VOXELIVIL 3 (ITIIII) | 1.25 | User def elem sel | no | Acquisition mode | cartesian |
| REC voxel MPS (mm) | 1.14 / 1.14 / | element selection | All | Fast Imaging mode | TFE |
| NEO VOXO: IIII O (IIIII) | 1.25 | connection | conn-A | 3D non-selective | no |
| Scan percentage (%) | 100 | Coil selection 2 | RX-Intf-2 | shot mode | multishot |
| TFE shots | 30 | element selection | All | TFE factor | |
| TFE dur. shot / acq (ms) | 1669.4 / 716.5 | Dual coil | yes | | 256 |
| Min. TI delay | 381.0045 | CLEAR | yes | 3D free factor | no default |
| Act. WFS (pix) / BW | 0.579 / 1750.7 | body tuned | yes | startup echoes | 0 |
| (Hz) | | FOV FH (mm) | 256 | +TFE followup echoes | user defined |
| Min. WFS (pix) / Max. | 0.575 / 1763.2 | AP (mm) | 256 | shot interval | |
| BW (Hz) | | RL (mm) | 172.5 | (ms) | 4500 |
| | 0.8993402 | Voxel size FH (mm) | 1.25 | profile order | linear |
| (W) | 22.04 | AP (mm) | 1.254902 | turbo direction | radial |
| SAR / head | < 33 % | RL (mm) | 1.254702 | CENTRA (spiral) | no |
| Whole body / level | 0.0 W/kg / | Recon voxel size (mm) | 1.142857 | Echoes | 1 |
| D1 | normal | Fold-over suppression | no | partial echo | no |
| B1 rms PNS / level // VUIIS : | 0.76 uT | Slice oversampling | no default | shifted echo | no |
| PNS / level // VUIIS : dortch : | 57 % / normal | RF select, FOS | | TE | shortest |
| Sound Pressure Level | 37.43258 | RE Select. FOS Reconstruction matrix | no 224 | Flip angle (deg) | 7 |
| (dB) | 37.43236 | l- | | TR | shortest |
| MOTION | J | SENSE (AD) | yes | Halfscan | no |
| | no | P reduction (AP) | 2 | Water-fat shift | minimum |
| Heart rate > 250 bpm | no | P os factor | 1 | Shim | auto |
| Respiratory | no | S reduction (RL) | 2 | mDIXON | no |
| compensation | 110 | k-t BLAST | no | Fat suppression | no |
| Navigator respiratory | no | Overcontiguous slices | no | Water suppression | no |
| comp | | Stacks | 1 | TFE prepulse | invert |
| Flow compensation | no | slices | 138 | slice selection | no |
| fMRI echo stabilisation | no | slice orientation | sagittal | delay | user defined |
| Motion smoothing | no | fold-over direction | AP | (ms) | 1300 |
| NSA | 1 | fat shift direction | F | PSIR | no |
| DYN/ANG | G | Chunks | 1 | +inv pulse type | +B1 opt (low |
| Angio / Contrast enh. | no | PlanAlign | no | | BW) |
| Quantitative flow | no | REST slabs | 0 | MTC | no |
| CENTRA | no | Catheter tracking | no | T2prep | no |
| Manual start | no | Interactive positioning | no | Research prepulse | no |
| +Abuse dynamic loop | no | Allow table movement | no | Diffusion mode | no |
| Dynamic study | no | OFFC/AN | G | Elastography mode | no |
| Arterial Spin labeling | no | Stacks | 1 | SAR mode | high |
| POST/PRO | | Stack Offc. AP | 11.5916 | B1 mode | default |
| Preparation phases | auto | (P=+mm) | | SAR Patient data | auto |
| Interactive F0 | | RL (L=+mm) | -3.434925 | PNS mode | low |
| SENSE ref. scan | no | FH (H=+mm) | 4.208135 | Gradient mode | full control |
| SmartPlan survey | | Ang. AP (deg) | 0 | max strength | 33 |
| B0 field map | no | RL (deg) | 0 | (mT/m) | |
| | no | FH (deg) | 0 | max slew rate | 166 |
| B1 field map MIP/MPR | no | l | | (T/m/s) | |
| | no | | | | |
| | no | | | | |
| Images | M, no, no, no | | | | |
| Images Autoview image | M, no, no, no | | | | |
| Images Autoview image Calculated images | M, no, no, no M no, no, no, no | | | | |
| Images Autoview image Calculated images Reference tissue | M, no, no, no M no, no, no, no Grey matter | | | | |
| Images Autoview image Calculated images Reference tissue Preset window contrast | M, no, no, no M no, no, no, no Grey matter soft | | | | |
| Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode | M, no, no, no M no, no, no, no Grey matter | | | | |
| Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode Save raw data | M, no, no, no M no, no, no, no Grey matter soft immediate no | | | | |
| Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode Save raw data Hardcopy protocol | M, no, no, no M no, no, no, no Grey matter soft immediate no no | | | | |
| Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode Save raw data Hardcopy protocol Ringing filtering | M, no, no, no M no, no, no, no Grey matter soft immediate no | | | | |
| Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode Save raw data Hardcopy protocol | M, no, no, no M no, no, no, no Grey matter soft immediate no no | | | | |

| Hospital (2) | |
|--|---------|
| Rel. signal level (%) 100 | |
| Act. TR/TE (ms) | |
| Dyn. scan time | |
| Time to k0 | |
| ACQ marrix M x P 160 x 148 | |
| ACQ voxel MPS (mm) 1.50 / 1.62 / 10.0 | |
| Connection Conn-A Consection Conn-A Connection Con | |
| Scan percentage (%) 92.77109 element selection Ali Mill | |
| Scan percentage (%) 92.77109 Dual coil yes Selective S | |
| Act. WFS (pix) / BW 8.715 / 116.3 CLEAR yes body tuned no FOV AP (mm) 240 FOV AP (mm) 240 FOV AP (mm) 240 FET (actor 7 For AP (mm) 1.5 FOV AP (mm) 1.5 FET (actor 7 For AP (actor 7 For AP (mm) 1.5 FET (actor 7 For AP | |
| Description | |
| BW in EPI freq. dir. (Hz) 1230.8 Min. WFS (pix) / Max 8.664 / 117.0 FOV AP (mm) 240 Min. TR/TE (ms) 54 / 5.9 RF avg power computed (W) SAR / head < 49 % RL (mm) 1.5 RR (mm) 1.5 | |
| Min. WFS (pix) / Max. 8.664 / 117.0 RL (mm) 240 direction RL (mm) 240 direction RL (mm) RL (mm) 165 Echoes 1 Mirch (mos) FH (mm) 165 Echoes 1 Mirch (mos) FH (mm) 1.5 Echoes 1 Mirch (mos) FH | |
| Min. TR/TE (ms) 54 / 5.9 RF avg power computed 1.351341 (W) Cover suppression 1.5 Recon voxel size (mm) 1.5 Recon voxel size (mm) 0.9375 Fold-over suppression no Filip angle (deg) 5 Reconstruction matrix 256 TR user defined Reconstruction matrix 256 TR user defined Respiratory Compensation no Cover outgoing size Navigator respiratory compensation NSA 1 FH (mm) 165 Echoes 1 Recon voxel size (mm) 0.9375 Fold-over suppression no Filip angle (deg) 5 Reconstruction matrix 256 TR user defined Resonstruction matrix 256 TR user defined Reconstruction matrix 256 TR user defined Reconstruction matrix 256 TR user defined Resonstruction matrix 256 TR user defined Reconstruction matrix 256 TR user defined Resonstruction matrix 256 TR | |
| MR MR MR MR MR MR MR MR | |
| RL (mm) 1.5 echo | |
| SAR / head | |
| Recon voxel size (mm) 0.9375 echo | |
| B1 rms 0.93 uT Slice oversampling default (ms) 7.2 | |
| Second S | |
| Reconstruction matrix 256 | |
| Sound Pressure Level (dB) MOTION Cardiac synchronization no Ferrorization Post factor 1 Water-fat shift minimum Not wolume Shim volume Not minimum | |
| Cardiac synchronization Preduction (RL) 2 Halfscan No | |
| MOTION P os factor 1 Water-fat shift minimum | |
| Cardiac synchronization no S reduction (FH) 2 Shim volume Heart rate > 250 bpm no Respiratory no compensation Navigator respiratory comp Flow compensation no Stacks 1 Fat suppression ProSet Slices 33 pulse type 1331 Slice orientation transverse slice orientation transverse fold-over direction RL fat shift direction L Chunks 1 Chunks 1 Chunks 1 PlanAlign no REST slabs 0 Hall or constant fold-over direction REST slabs 1 Catheter tracking no Hall or constant fold-over direction REST slabs 1 Chunks | |
| Heart rate > 250 bpm no Respiratory no Compensation Navigator respiratory comp no Stacks 1 Fat suppression ProSet Stacks 1 Stacks 1 Suppression ProSet Stacks 1 Stacks 1 Suppression ProSet | |
| Respiratory compensation Navigator respiratory comp Flow compensation NAI echo stabilisation NA DYN/ANG Anglo / Contrast enh. Quantitative flow Anaual start +Abuse dynamic loop Dynamic study dyn scans FoV time mode dummy scans immediate subtraction NA Overcontiguous slices yes MDIXON NA Fat suppression ProSet Silces 33 Silce orientation Itransverse Indid over direction RL MTC HOURT MAIC HOURT Chunks 1 Ch | |
| Stacks 1 | |
| Silce orientation Transverse Water Suppression Silce orientation Silce orientation Transverse Water Suppression Silce orientation Silce orientation Transverse Water Suppression Silce orientation | |
| Silice orientation Transverse Suppression | |
| FMRI echo stabilisation no NSA | |
| The image of the | OT |
| Angio / Contrast enh. no PlanAlign no REST slabs 0 +Amp 3 | :51 |
| Angio / Contrast enh. no Quantitative flow no Guantitative flow no Amanual start no Hahual start no Hahual start no Dynamic study individual dyn scans 64 recon multiplier 1 dyn scan times shortest FOV time mode dummy scans immediate subtraction FIGURE FLOW Stack Stack Stack Stack Offic. AP CP + mm) RL (L=+mm) -3.434925 FH (H=+mm) 35.457352 Ang. AP (deg) 0 Research no PlanAlign no HB1 Mode constant +Amp 3 +B1 Units max amp. (uT) +Offset file Mode FURTHER STACK G:/patch/rf_off File +Offset Vinits +Offset File +Offset Vinits +Offset Popm +Offset Stack G:/patch/rf_off File +Offset Popm +Offset Stack G:/patch/rf_off File +Offset Popm +Offset Stack G:/patch/rf_off File +Offset Popm +Offset | |
| Manual start no Manual start no HADUSE dynamic loop no Dynamic study individual dyn scans FOV time mode dummy scans 0 immediate subtraction Manual start no Interactive positioning no Allow table movement no OFFC/ANG Stack offc. AP (P=+mm) RL (L=+mm) 43.434925 FH (H=+mm) 35.45387 Ang. AP (deg) 0 Research no Tenulus max amp. (uT) +Pulses/TR 1 +Offset file Mode 1-FOFSet File +Offset Units FILE +Offset Pl +Offset Units PM -OFFEC/ANG Stack offc. AP (P=+mm) -3.434925 FH (H=+mm) 35.45387 Ang. AP (deg) 0 Research no Research no Tenulus max amp. (uT) +Pulses/TR 1 +Offset File +Offset Units -OFFEC/ANG Stack offc. AP (P=+mm) -3.434925 FH (H=+mm) 35.45387 Ang. AP (deg) 0 Research no Tenulus max amp. (uT) | |
| Manual start no | |
| Dynamic study individual dyn scans 64 recon multiplier 1 dyn scan times shortest FOV time mode dummy scans immediate subtraction Allow table movement no +Offset file Mode Stack S 1 +Offset File +Offset G:/patch/rf_off File +Offset Dunits Photoset File +Offset Dunits Photoset Photoset File +Offset Dunits Photoset Photoset File +Offset Dunits Photoset Photose | |
| Dynamic study Individual OFFC/ANG Mode +Offset file | |
| Stack 1 Stack Stac | |
| Stack Offc. AP | oto tut |
| (P=+mm) -3.434925 +Offset Units +Offset +Off | eis.ixi |
| dummy scans 0 RL (L=+mm) -3.434925 Units immediate subtraction no Ang. AP (deg) 0 Research no repulse | |
| immediate subtraction PI (H=+mm) 35.45897 +RF Shape gauss_cest_mt Ang. AP (deg) 0 Research no repulse | |
| subtraction Ang. AP (deg) 0 Research no | |
| | |
| Tast next scan I no | |
| synch. ext. device no Shim Sizo AD (mm) 190 7405 Flastography no | |
| MIC no Pl (mm) 134 1432 mode | |
| dyn stabilization no FH (mm) 114 2944 SAR mode high | |
| offic. AP (P=+mm) -1.87074 B1 mode default | |
| RI (I = +mm) -3.660397 SAR Patient auto | |
| Arterial Spin labeling no FH (H=+mm) 33.86363 | |
| POST/PROC Ang. AP (deg) 2.26406 First mode default | |
| Preparation phases sameprep RL (deg) -6.580403 SofTone mode no | |
| Interactive FO no FH (deg) 4.637685 | |
| SENSE ref. scan no | |
| SmartPlan survey no | |
| B0 field map no | |
| B1 field map no | |
| MIP/MPR no Images M, no, no, no | |
| Autoview image M | |
| Calculated images no, no, no, no | |
| Reference tissue Grey matter | |
| EPI 2D phase correction no | |
| Preset window contrast soft | |
| Reconstruction mode real time | |
| reuse memory no | |
| Save raw data no | |
| Hardcopy protocol no | |
| Ringing filtering rectangular | |
| Geometry correction default | |
| Elliptical k-space shutter default | |

| ☐ Hospital (2) ☐ 201 | 41021 CEST fMRI | (13) 52:37.4 CEST_ | Reddy_GluCEST | 11:37.0 | |
|---|--|---|---------------------|--------------------------|---------------------|
| INFO PA | | GEOMETR | | CONTRA | ST |
| Total scan duration | 11:37.0 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | 2D |
| Act. TR/TE (ms) | 5.6 / 2.7 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| Dyn. scan time | 00:13.909 | Xmit Coil selection | MTX-Volume- | + ZOOM | no |
| Time to k0 | 00:08.6 | User def elem sel | T/R | Contrast enhancement | T1 |
| ACQ matrix M x P | 128 x 126 | element selection | no All | Acquisition mode | cartesian |
| ACQ voxel MPS (mm) | 1.88 / 1.90 / 10.0 | connection | conn-A | Fast Imaging mode | TFE |
| REC voxel MPS (mm) | 0.94 / 0.94 / | Coil selection 2 | RX-Intf-2 | shot mode TFE factor | multishot 3 |
| The voxor in a (inin) | 10.0 | element selection | All | startup echoes | default |
| Scan percentage (%) | 98.4375 | Dual coil | yes | +TFE followup echoes | 0 |
| TFE shots | 42 | CLEAR | yes | shot interval | shortest |
| TFE dur. shot / acq (ms) | | body tuned | no | profile order | low_high |
| TFE shot interval (ms) | 331.1945 | FOV AP (mm) | 240 | Echoes | 1 |
| Act. WFS (pix) / BW (Hz) | 0.700 / 1446.8 | RL (mm) | 240 | partial echo | no |
| Min. WFS (pix) / Max. | 0.649 / 1562.5 | Voxel size AP (mm) | 1.87 | shifted echo | no |
| BW (Hz) | | RL (mm) | 1.87 | TE | user defined |
| Min. TR/TE (ms) | 4.0 / 1.36 | Slice thickness (mm) | 10 0.9375 | (ms) | 2.7 |
| RF avg power computed | 2.058919 | Recon voxel size (mm) Fold-over suppression | no | Flip angle (deg) | 10 |
| (W) | . 75 0/ | Reconstruction matrix | 256 | TR (ms) | user defined 5.6 |
| SAR / head | < 75 % | SENSE | no | (ms) Halfscan | 5.6 no |
| Whole body / level | < 0.1 W/kg / normal | k-t BLAST | no | Water-fat shift | minimum |
| B1 rms | 1.15 uT | Slice orientation | transverse | Shim | volume |
| PNS / level // VUIIS : | 31 % / normal | Fold-over direction | RL | ShimAlign | no |
| dortch : | | Fat shift direction | P | mDIXON | no |
| Sound Pressure Level | 17.13024 | PlanAlign | no | Fat suppression | no |
| (dB) | <u> </u> | REST slabs | 0 | Water suppression | no |
| MOTIO | | Catheter tracking | no | TFE prepulse | no |
| Cardiac synchronization Heart rate > 250 bpm | no no | Interactive positioning | no | MTC | + pulsed |
| Respiratory | no | Allow table movement | no | | qMT/CEST |
| compensation | 110 | OFFC/AN | | +Duration (ms) | 10 |
| Navigator respiratory | no | Slice Offc. AP (P=+mm) | -4.247866 | +B1 Mode +Amp | constant 4.25 |
| comp | | RL (L=+mm) | -3.434925 | +Amp +B1 Units | max amp. (uT) |
| Flow compensation | no | FH (H=+mm) | 35.88707 | +Pulses/TR | 10 |
| fMRI echo stabilisation | no | Ang. AP (deg) | 0 | + Duty Cycle | 0.9 |
| Motion smoothing | no | RL (deg) | -5.457352 | +Offset Mode | baseline+range |
| NSA | 1 | FH (deg) | 0 | +Min Offset | -5 |
| DYN/AN Angio / Contrast enh. | no | Shim Size AP (mm) | 180.7495 | +Max Offset | 5 |
| Quantitative flow | no | RL (mm) | 134.1432 | +Offset Units | ppm |
| Manual start | no | FH (mm) | 106.6376 | +RF Shape | gauss_cest_mt |
| +Abuse dynamic loop | no | Offc. AP (P=+mm) | -1.87074 | + Interpulse | no |
| Dynamic study | individual | RL (L=+mm) | -3.660397 | Spoiling | |
| dyn scans | 50 | FH (H=+mm) Ang. AP (deg) | 33.86363 2.26406 | T2prep Research prepulse | no no |
| recon multiplier | 1 | RL (deg) | -6.580403 | Diffusion mode | no |
| dyn scan times | shortest | FH (deg) | 4.637685 | Elastography mode | no |
| FOV time mode | default | (dog) | 1.007000 | SAR mode | moderate |
| dummy scans | 0 | | | B1 mode | default |
| immediate subtraction | no | | | SAR Patient data | auto |
| fast next scan | no | 1 | | PNS mode | low |
| synch. ext. device | no | † | | Gradient mode | default |
| MTC MTC | no | † | | SofTone mode | no |
| dyn stabilization | no | † | | | |
| prospect. motion | no | | | | |
| corr. | | | | | |
| Keyhole | no | | | | |
| Arterial Spin labeling | no | | | | |
| | UC | ļ | | | |
| Propagation phases | camonron | | | | |
| Preparation phases | sameprep | | | | |
| Preparation phases Interactive F0 | no | | | | |
| Preparation phases Interactive F0 SENSE ref. scan | no no | | | | |
| Preparation phases Interactive F0 | no | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey | no no no | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map | no no no | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map | no no no no no | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image | no no no no no no | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Ilm/MPR Ilm/MPS Autoview image Calculated images | no no no no no M, no, no, no M | | | | |
| Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | no no no no no no M, no, no, no M ro, no, no, no Grey matter | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast | no no no no no no no M, no, no, no M Grey matter | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode | no on M, no, no, no, no M on, no, no, no Grey matter soft real time | | | | |
| Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode reuse memory | no no no no no no no M, no, no, no M no, no, no, no Grey matter soft real time no | | | | |
| Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode reuse memory Save raw data | no no no no no no mo M, no, no, no M no, no, no, no Grey matter soft real time no no | | | | |
| Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode reuse memory Save raw data Hardcopy protocol | no no no no no no no M, no, no, no, no M or, no, no, no, no Grey matter soft real time no no no | | | | |
| Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode reuse memory Save raw data | no no no no no no mo M, no, no, no M no, no, no, no Grey matter soft real time no no | | | | |

| Modern 1962 | Hospital (2) 201 | 41021 CEST fMRI | (13) 52:37 4 FR1 Re | ddy multi∆nale∩ | 11.42 0 | |
|--|---|-------------------------|---------------------------------------|-----------------|--|--------------|
| Total scan duration 0.14.2 0 Muldius minum no | 1 | | | | | ет |
| Ref. Signal level (%) | | | | | | |
| Act. TE (mm) | | | | | | |
| Act. Tet (ren) 12 2 | | | | | | |
| Dyn. scan filme | <u> </u> | | | | | |
| ACO mortifix M. F. (128 x 120) User def elem sell on olements selection of elements selection of elements selection of the elements of elements selection of the elements of the elements selection of the elements of the eleme | | | | | | - |
| 10.0 | | | User def elem sel | no | | no |
| Connection Con | | | element selection | All | | cartesian |
| San percentage (%) 93.75 Dual of II West (%) 93.75 Dual of II West (%) 93.75 Dual of II West (%) 10 CLEAR No | ` ′ | 10.0 | connection | conn-A | | TSE |
| Scan percentage (%) 3.75 She percentage (%) 10 CICAP (mm) 240 Final (mm) 10 CicAP (mm) 240 Final (mm) 2 | REC voxel MPS (mm) | | Coil selection 2 | RX-Intf-2 | shot mode | multishot |
| Packages | | | l | All | TSE factor | 15 |
| Min. size gap (mm) 10 10 10 10 10 10 10 1 | | | | - | startup echoes | 0 |
| Wis Sigh / BW (Hz) 11.986 / \$10.1 Rf. (mm) 240 DRIVE NoMyl. TEEF / Tecquiv (ms) 12.2 / 12 12.9 / 130 FH. (mm) 10 utrashed no Min. TR (ms) 589 RF. (mm) 1.87 country country country RR. (mm) 1.87 country country country no country country no country no country no month no no no country no month no no month no no month no month no no month no no <td< td=""><td></td><td></td><td>I</td><td></td><td></td><td>0</td></td<> | | | I | | | 0 |
| TSE est / Shot (ms) 12 0 / 180 FH (mm) 10 Ultrahort 10 Ultra | | | | | | low_high |
| Teef Teequiv (ms) 12 / 12 12 18 18 18 18 18 18 | | | | | | no |
| Min. TR (ms) S89 | | | (| | | |
| Rice https://doi.org/10.1006 | | | | - | | no |
| No. | | | | | | 1 |
| SAR / head 1 0 % Small FOV imaging no TE user defined Whole body / level 0.0 4 l/ u T Phys / level // Yull S: 0.4 l/ u T Phys / level // Yull S: 30 % / normal Reconstruction matrix 256 Flip angle (deg) 0.0 Sound Pressure Level (dis) 16.02466 Stacks no Refound (ms) 10 ms 4 Ref Pulse Dur (msc) 4 Ref Pulse Dur (msc) 4 Stacks 1 ms (ms) 6000 4 Ref MIX Dephase 5 Stacks 1 ms Water-fat shift user defined (glee) (slee) 2 2 4 Ref MIX Dephase 5 Stacks 1 ms Water-fat shift user defined (glee) (slee) 2 2 4 Ref MIX Dephase 5 Stacks 1 ms Stacks 1 ms Water-fat shift user defined (glee) (slee) (slee) 2 2 4 Ref MIX Dephase 5 Stacks 1 ms Stack fill drection R R. ms ms Mill pressor No Shift drection R R. mp ps Shift drection R R. ms ps ps | | | | | + | |
| More Do Wins normal No | SAR / head | < 10 % | | | | |
| BT rms | Whole body / level | | | | | |
| SFMS Tower // VUIIS 30 % normal dortch | D1 | | | | ` ' | |
| Manual start | | ***** | | | | |
| Stacks 1 | | 30 % / normal | \ <u> </u> | | | |
| Subsect Subs | | 16.02466 | | | + | |
| Silces 1 | | 10.02400 | | | | no |
| Ref Pulse Dur [msec] | · / | | | | Water-fat shift | user defined |
| Shim Volume Vol | | | slice gap | user defined | (pixels) | 2 |
| Activation Cardiac synchronization No P Fat suppression No P | | | | 0 | Shim | volume |
| MOTION Cardiac synchronization no Minimum number of packages Sice scan order Sice scan order Alphanaging no Minimum number of packages Grad. rev. offres. supp. no Navigator respiratory no Cardiac spacing Minimum number of planAlign no Minimum number of packages Minimum nu | +Ref MIN Dephase | | slice orientation | transverse | ShimAlign | no |
| Cardiac synchronization no Heart rate > 250 bpm no Silice scan order default Respiratory compensation no Mayagator respiratory comp Navigator respiratory no Catheter tracking no Catheter tracking no Temporal silice spacing Motion smoothing no Motion smoothing no NSA 1 Stack Offc. AP -4.247866 (P=+mm) -3.434925 FH (H=+mm) 35.88707 Ang. AP (deg) 0 My scans 2 PH (deg) 0 My scans 0 FH (mm) 134.1432 My scans 0 FH (deg) 0 0 0 0 0 0 0 0 0 | +Crusher b value | | fold-over direction | RL | | no |
| Deckages Deckages Silce scan order Deckages D | MOTIO | N | \ <u> </u> | · | | no |
| Respiratory compensation no Respiratory compensation no Research preparation no Research prepulse no Catheter tracking no Interactive positioning no Research prepulse no Catheter tracking no Interactive positioning no Research prepulse no DYN/ANG Allow table movement no Diffusion mode no Diffusion m | | no | | 1 | | |
| PlanAlign no | | | | dofault | | |
| Navigator respiratory comp Comp Comp Comp Comp Comp Comp Comp C | | no | l—————— | | | |
| Catheter tracking no Research prepulse no Tomporal slice spacing default Molton smoothing no NSA | | | | | | |
| Flow compensation no | | no | | | | |
| Allow table movement no Diffusion mode no | | no | | - | | |
| Motion smoothing no NSA 1 Stacks 1 Stack Offic. AP -4.247866 Manual start no +Abuse dynamic loop diminish minish | | | | | | |
| NSA | | | - | | | |
| Namual start | | 1 | | | | |
| Manual start | DYN/AN | G | | | | - |
| +Abuse dynamic loop diminish Dynamic study individual Agn Scans 2 PH (H=+mm) 35.88707 Dynamic study individual Agn AP (deg) 0 PNS mode Nigh Gradient mode default Agn Scans 1 PH (H=+mm) 35.88707 Fecon multiplier 1 FH (H=+mm) 35.88707 FH (deg) 0 PHS mode Mefault Agn PHS mode Mefault Agn PHS mode Mefault Agn PHS Mega | Manual start | no | (P=+mm) | | | |
| Dynamic study individual dyn scans 2 RL (deg) -5.457352 Second tribulation of the state of the s | +Abuse dynamic loop | | | | | auto |
| dyn scans 2 RL (deg) -5.457352 SofTone mode default FH (deg) 0 SofTone mode no SofTone mode sofTone | | | · · · · · · · · · · · · · · · · · · · | | PNS mode | high |
| recon multiplier 1 dyn scan times shortest FOV time mode default dummy scans 0 FOV time mode default no subtraction | | | | | Gradient mode | default |
| dyn scan times shortest FOV time mode default dummy scans 0 immediate subtraction fast next scan no synch. ext. device dyn stabilization no prospect. motion corr. Keyhole no Arterial Spin labeling no POST/PROC Preparation phases sameprep Interactive FO no SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Autoview image M RL (mm) 134.1432 FH (mm) 106.6376 FH (H=+mm) -3.660397 FH (H=+mm) 33.86363 Ang. AP (deg) 2.26406 RL (deg) -6.580403 FH (deg) 4.637685 FH (deg) 4.637685 | | | . 0, | | SofTone mode | no |
| FOV time mode default dummy scans 0 | | | | | | |
| dummy scans 0 FH (mm) 106.6376 immediate no subtraction | | | | | | |
| immediate subtraction | | | | | | |
| subtraction fast next scan no | | | | | | |
| fast next scan no SHL (H=+mm) 33.86363 synch. ext. device no Ang. AP (deg) 2.26406 dyn stabilization no prospect. motion corr. Keyhole no Arterial Spin labeling no POST/PROC Preparation phases sameprep Interactive FO no SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | 1 | | | | |
| synch. ext. device no Ang. AP (deg) 2.26406 dyn stabilization no prospect. motion corr. Keyhole no Arterial Spin labeling no POST/PROC Preparation phases sameprep Interactive FO no SENSE ref. scan no SmartPlan survey no B0 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | no | | | | |
| dyn stabilization no RL (deg) -6.580403 prospect. motion no RL (deg) -6.580403 FH (deg) - | | | | | | |
| prospect. motion corr. Keyhole no Arterial Spin labeling no POST/PROC Preparation phases sameprep Interactive FO no SENSE ref. scan no SmartPlan survey no B0 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Reference tissue Grey matter Preset window contrast Soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | no | | | | |
| corr. Keyhole no Arterial Spin labeling no POST/PROC Preparation phases sameprep Interactive FO no SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no Images M, no, no, no Autoview image M Calculated images no, | | no | | | | |
| Arterial Spin labeling no POST/PROC Preparation phases sameprep Interactive F0 no SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Ringing filtering rectangular | | | (| | | |
| POST/PROC Preparation phases sameprep Interactive FO no SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no Images M, no, no, no Autoview image M Calculated images no, | | | | | | |
| Preparation phases sameprep Interactive F0 no SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Reference tissue Grey matter Preset window contrast Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | | | | |
| Interactive FO no SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Ringing filtering rectangular | | | | | | |
| SENSE ref. scan no SmartPlan survey no B0 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Ringing filtering rectangular | | | | | | |
| SmartPlan survey no B0 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | | | | |
| B0 field map no B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | } | | | |
| B1 field map no MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | } | | | |
| MIP/MPR no Images M, no, no, no Autoview image M Calculated images no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | | | | |
| Images M, no, no, no Autoview image M Calculated images no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | } | | | |
| Autoview image M Calculated images no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | | | | |
| Calculated images no, no, no, no Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | † | | | |
| Reference tissue Grey matter Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | İ | | | |
| Preset window contrast soft Reconstruction mode real time reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | | | | |
| reuse memory no Save raw data no Hardcopy protocol no Ringing filtering rectangular | | | 1 | | | |
| Save raw data no Hardcopy protocol no Ringing filtering rectangular | | real time | | | | |
| Hardcopy protocol no Ringing filtering rectangular | | no | | | | |
| Ringing filtering rectangular | reuse memory | 110 | | | | |
| | Save raw data | no | | | | |
| Geometry correction default | Save raw data Hardcopy protocol | no no | | | | |
| | Save raw data Hardcopy protocol Ringing filtering | no no rectangular | | | | |

| ☐ Hospital (2) ☐ 201 | 41021 CEST fMRI | (13) 52:37.4 B0_Re | ddy_mulitecho 00 | 0:03.9 | |
|-------------------------|------------------------|-------------------------------|------------------|----------------------|-----------|
| INFO PAG | GE | GEOMETR | !Y | CONTRA | ST |
| Total scan duration | 00:03.9 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | 2D |
| Act. TR/TE1/delta TE | 53 / 3.4 / 3.9 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| (ms) | | Xmit Coil selection | MTX-Volume- | + ZOOM | no |
| ACQ matrix M x P | 128 x 128 | | T/R | Contrast enhancement | no |
| ACQ voxel MPS (mm) | 1.88 / 1.88 / 10.0 | User def elem sel | no | Acquisition mode | cartesian |
| REC voxel MPS (mm) | 0.94 / 0.94 / | element selection | All | Fast Imaging mode | none |
| REC VOXELIVIES (ITIIII) | 10.0 | connection | conn-A | Echoes | 4 |
| Scan percentage (%) | 100 | Coil selection 2 | RX-Intf-2 | partial echo | no |
| Act. WFS (pix) / BW | 2.824 / 358.8 | element selection | All | shifted echo | no |
| (Hz) | | Dual coil CLEAR | yes | TE first | shortest |
| Min. WFS (pix) / Max. | 0.649 / 1562.5 | <u> </u> | yes | echospacing | shortest |
| BW (Hz) | | body tuned | no | flyback | yes |
| RF avg power computed | 2.058916 | FOV AP (mm) | 240 | Flip angle (deg) | 65 |
| (W) | 0. | RL (mm) | 240 | TR | shortest |
| SAR / head | < 75 % | Voxel size AP (mm) | 1.87 | Halfscan | no |
| Whole body / level | < 0.1 W/kg / normal | RL (mm) | 1.875 | Water-fat shift | maximum |
| B1 rms | 1.15 uT | Slice thickness (mm) | 10 | Shim | volume |
| PNS / level // VUIIS : | 46 % / normal | Recon voxel size (mm) | 0.9375 | ShimAlign | no |
| dortch : | 40 % / 110111141 | Fold-over suppression | no | mDIXON | no |
| Sound Pressure Level | 21.87152 | Reconstruction matrix | 256 | Fat suppression | no |
| (dB) | 21.07.102 | SENSE | yes | Water suppression | no |
| MOTIO | N | P reduction (RL) | 2 | MTC | no |
| Cardiac synchronization | no | P os factor | 1 | Research prepulse | no |
| Heart rate > 250 bpm | no | k-t BLAST | no . | Diffusion mode | no |
| Respiratory | no | Slice orientation | transverse | Elastography mode | no |
| compensation | | Fold-over direction | RL | SAR mode | moderate |
| Navigator respiratory | no | Fat shift direction PlanAlign | P no | B1 mode | default |
| comp | | REST slabs | 0 | SAR Patient data | auto |
| Flow compensation | yes | Catheter tracking | no | PNS mode | low |
| fMRI echo stabilisation | no | Interactive positioning | no | Gradient mode | default |
| NSA | 1 | Allow table movement | no | SofTone mode | no |
| DYN/AN | | OFFC/AN | | | |
| Angio / Contrast enh. | no | Slice Offc. AP | -4.247866 | | |
| Quantitative flow | no | (P=+mm) | -4.24/000 | | |
| Manual start | no | RL (L=+mm) | -3.434925 | | |
| +Abuse dynamic loop | no | FH (H=+mm) | 35.88707 | | |
| Dynamic study | no | Ang. AP (deg) | 0 | | |
| Arterial Spin labeling | no | RL (deg) | -5.457352 | | |
| POST/PR | | FH (deg) | 0 | | |
| Preparation phases | sameprep | Shim Size AP (mm) | 180.7495 | | |
| Interactive F0 | no | RL (mm) | 134.1432 | | |
| SENSE ref. scan | no | FH (mm) | 106.6376 | | |
| SmartPlan survey | no | Offc. AP (P=+mm) | -1.87074 | | |
| B0 field map | no | RL (L=+mm) | -3.660397 | | |
| B1 field map | no | FH (H=+mm) | 33.86363 | | |
| MIP/MPR | no | Ang. AP (deg) | 2.26406 | | |
| Images | M, R, I, no | RL (deg) | -6.580403 | | |
| Autoview image | M no no no | FH (deg) | 4.637685 | [| |
| Calculated images | no, no, no, no | | | • | |
| Reference tissue | Grey matter | | | | |
| Preset window contrast | soft | | | | |
| Reconstruction mode | real time | | | | |
| Save raw data | yes | | | | |
| Hardcopy protocol | no | | | | |
| Ringing filtering | rectangular | | | | |
| Geometry correction | default | <u> </u> | | | |

| | | (13) 52:37.4 qMT H | | 0001=== | rT. |
|--|--|-------------------------|-------------|---|---|
| INFO PAG | | GEOMETR | | CONTRA | <u> </u> |
| Total scan duration | 02:38.7 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | 3D |
| Act. TR/TE (ms) | 4.2 / 2.2 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| Dyn. scan time | 00:11.251 | Xmit Coil selection | MTX-Volume- | + ZOOM | no |
| Time to k0 | 00:07.0 | ļ | T/R | Contrast enhancement | T1 |
| ACQ matrix M x P | 212 x 210 | User def elem sel | no | Acquisition mode | cartesian |
| ACQ voxel MPS (mm) | 1.00 / 1.01 / | element selection | All | Fast Imaging mode | TFE |
| | 2.00 | connection | conn-A | 3D non-selective | no |
| REC voxel MPS (mm) | 0.95 / 0.95 / | Coil selection 2 | RX-Intf-2 | shot mode | multishot |
| . (01) | 2.00 | element selection | All | TFE factor | 54 |
| Scan percentage (%) | 99.08257 | Dual coil | yes | 3D free factor | no |
| TFE shots | 12 | CLEAR | yes | startup echoes | user defined |
| TFE dur. shot / acq (ms) | | body tuned | no | (number) | 0 |
| TFE shot interval (ms) | 3353.226 | FOV AP (mm) | 212 | +TFE followup echoes | 0 |
| Act. WFS (pix) / BW (Hz) | 1.142 / 887.5 | RL (mm) | 212 | shot interval | shortest |
| Min. WFS (pix) / Max. | 1.133 / 894.1 | FH (mm) | 10 | profile order | low_high |
| BW (Hz) | 1.133 / 694.1 | Voxel size AP (mm) | 1 | turbo direction | Υ |
| RF avg power computed | 1 347757 | RL (mm) | 1 | Echoes | 1 |
| (W) | 1.547757 | FH (mm) | 2 | partial echo | no |
| SAR / head | < 49 % | Recon voxel size (mm) | 1 | shifted echo | no |
| Whole body / level | < 0.1 W/kg / | Fold-over suppression | no | TE | shortest |
| | normal | Slice oversampling | default | Flip angle (deg) | 15 |
| B1 rms | 0.93 uT | RF select. FOS | no | TR | shortest |
| PNS / level // VUIIS : | 50 % / normal | Reconstruction matrix | 224 | Halfscan | no |
| dortch : | | SENSE | yes | Water-fat shift | minimum |
| Sound Pressure Level | 29.07384 | P reduction (RL) | 2 | Shim | PB-volume |
| (dB) | | P os factor | 1 | ShimAlign | no |
| MOTIO | N | S reduction (FH) | 1 | mDIXON | no |
| Cardiac synchronization | no | k-t BLAST | no | Fat suppression | no |
| Heart rate > 250 bpm | no | Overcontiguous slices | no | Water suppression | no |
| Respiratory | no | Stacks | 1 | TFE prepulse | no |
| compensation | | slices | 5 | | |
| Navigator respiratory | no | slice orientation | transverse | MTC | +SIR |
| comp | | fold-over direction | RL | +pulse dur (ms) | 5.5 |
| Flow compensation | no | fat shift direction | P | +pulse shape | BRASORF |
| fMRI echo stabilisation | no | Chunks | 1 | +offset (Hz) | 0 |
| Motion smoothing | no | PlanAlign | no | +spoil amp (mT/m) | 20 |
| NSA | 1 | REST slabs | 0 | +spoil dur (ms) | 2 |
| DYN/AN | G | Catheter tracking | no | +TFE Saturation | |
| Angio / Contrast enh. | no | | no | +td Mode | yes constant |
| Quantitative flow | no | Interactive positioning | | <u> </u> | 2500 |
| CENTRA | no | Allow table movement | no | +td (ms) | |
| Manual start | no | OFFC/AN | | +ti (ms) | 6, 10, 16, 26, 42, 68, 110, |
| +Abuse dynamic loop | no | Stacks | 1 | | 178, 288, 468 |
| Dynamic study | individual | Stack Offc. AP | -4.247866 | | 760, 1233, |
| dyn scans | 14 | (P=+mm) | 2.424025 | | 2000, 8000, |
| recon multiplier | 1 | RL (L=+mm) | -3.434925 | | 2000, 10000, |
| dyn scan times | shortest | FH (H=+mm) | 35.88707 | | 1000, 1000, 1000, 1000, |
| FOV time mode | default | Ang. AP (deg) | 0 | | 1000, 1000, |
| dummy scans | 0 | RL (deg) | -5.457352 | | 1000, 1000, |
| immediate | no | FH (deg) | 0 | | 1000, 1000, |
| subtraction | | Shim Size AP (mm) | 180.7495 | | 1000, 1000, |
| fast next scan | no | RL (mm) | 134.1432 | | 1000, 1000, |
| synch. ext. device | no | FH (mm) | 106.6376 | | 1000, 1000, 1000, 1000, |
| MTC | no | Offc. AP (P=+mm) | | | 1000, 1000, |
| dyn stabilization | no | RL (L=+mm) | -3.660397 | | 1000, 1000, |
| prospect. motion | no | FH (H=+mm) | 33.86363 | | 1000, 1000, |
| corr. | - | Ang. AP (deg) | 2.26406 | | 1000, 1000, 1000, 1000, |
| Keyhole | no | RL (deg) | -6.580403 | | 1000, 1000, |
| Arterial Spin labeling | no | FH (deg) | 4.637685 | | 1000, 1000, |
| | | | | | 1000, 1000, |
| POST/PRO | oc | | | i . | 1000, 1000, |
| POST/PRO | | | | | 4000 |
| POST/PRO Preparation phases | auto | | | | 1000, 1000, |
| POST/PRO Preparation phases Interactive F0 | auto no | | | | 1000, 1000, |
| POST/PRO Preparation phases Interactive F0 SENSE ref. scan | auto no no | | | | |
| POST/PRO Preparation phases Interactive FO SENSE ref. scan SmartPlan survey | no no no | | | | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, |
| POST/PRO Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map | auto no no no no | | | | 1000, 1000, 1000, 1000, 1000, 1000, |
| POST/PRO Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map | auto no no no no no no | | | T2prep | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, |
| POST/PRO Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR | auto no | | | T2prep Research prepulse | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 |
| POST/PRO Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images | auto no no no no no no no M, R, I, no | | | | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 |
| POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image | auto no no no no no no no M, R, I, no M | | | Research prepulse | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no |
| POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images | auto no no no no no no M, R, I, no M no, no, no, no, no | | | Research prepulse Diffusion mode Elastography mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no |
| POST/PRd Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | auto no no no no no no no M, R, I, no M | | | Research prepulse Diffusion mode Elastography mode SAR mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no |
| POST/PRd Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | auto no no no no no no M, R, I, no M no, no, no, no, no | | | Research prepulse Diffusion mode Elastography mode SAR mode B1 mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no no high default |
| POST/PRe Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast | auto no no no no no no no no M, R, I, no M no, no, no, no, no Grey matter | | | Research prepulse Diffusion mode Elastography mode SAR mode B1 mode SAR Patient data | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no no high default auto |
| | auto no no no no no no no M, R, I, no M on, no, no, no, no Grey matter | | | Research prepulse Diffusion mode Elastography mode SAR mode B1 mode SAR Patient data PNS mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no no high default auto |
| POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode | auto no no no no no no no no M, R, I, no M on, no, no, no, no, no Grey matter soft real time | | | Research prepulse Diffusion mode Elastography mode SAR mode B1 mode SAR Patient data PNS mode Gradient mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no high default auto high |
| POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode reuse memory Save raw data | auto no M, R, I, no M no, no, no, no, no Grey matter soft real time no | | | Research prepulse Diffusion mode Elastography mode SAR mode B1 mode SAR Patient data PNS mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no high default auto |
| POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode reuse memory Save raw data Hardcopy protocol | auto no no no no no no no no no M, R, I, no M no, no, no, no, no Grey matter soft real time no no no | | | Research prepulse Diffusion mode Elastography mode SAR mode B1 mode SAR Patient data PNS mode Gradient mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no high default auto high |
| POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue Preset window contrast Reconstruction mode reuse memory | auto no no no no no no no no M, R, I, no M no, no, no, no, no Grey matter soft no no no | | | Research prepulse Diffusion mode Elastography mode SAR mode B1 mode SAR Patient data PNS mode Gradient mode | 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000, 1000 no no no high default auto high |

| Hospital (2) 201 | A1021 CEST MADE | (13) 52:37 / TEMPI | | 08:34 0 | |
|---|--|---------------------------|--|--------------------------|--------------|
| INFO PAG | | ii ee | _RESTINGSTATE | 11 | ет |
| Total scan duration | 08:34.0 | GEOMETE Multi-transmit | | Scan type | Imaging |
| | | | no H1 | Scan type | |
| Rel. signal level (%) | 100 | Nucleus | | Scan mode | MS |
| Act. TR/TE (ms) | 2000 / 25 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| Dyn. scan time | 00:02.000 | Xmit Coil selection | MTX-Volume- T/R | + ZOOM | no |
| Time to k0 | 00:15.0 | User def elem sel | no | Contrast enhancement | no |
| ACQ matrix M x P | 96 x 95 | l | | Acquisition mode | cartesian |
| ACQ voxel MPS (mm) | 2.50 / 2.50 / | element selection | All | Fast Imaging mode | EPI |
| | 2.50 | connection | conn-A | shot mode | single-shot |
| REC voxel MPS (mm) | 2.50 / 2.50 / | Coil selection 2 | RX-Intf-2 | Echoes | 1 |
| . (01) | 2.50 | element selection | All | partial echo | no |
| Scan percentage (%) | 100 | Dual coil | yes | shifted echo | no |
| Packages | 1 | CLEAR | yes | TE | user defined |
| Min. slice gap (mm) | 0 | body tuned | yes | (ms) | 25 |
| EPI factor | 37 | FOV RL (mm) | 240 | Flip angle (deg) | 63 |
| Act. WFS (pix) / BW | 21.825 / 46.4 | AP (mm) | 240 | TR | user defined |
| (Hz) | | FH (mm) | 115 | (ms) | 2000 |
| BW in EPI freq. dir. (Hz) | 2878.9 | Voxel size RL (mm) | 2.5 | Halfscan | no |
| Min. WFS (pix) / Max. | 21.770 / 46.5 | AP (mm) | 2.5 | Water-fat shift | - |
| BW (Hz) | | Slice thickness (mm) | 2.5 | 7 | minimum |
| Min. TR/TE (ms) | 1999 / 12 | Recon voxel size (mm) | 2.5 | Shim | auto |
| RF avg power computed | 1.371604 | l | | mDIXON | no |
| (W) | | Fold-over suppression | no | Fat suppression | no |
| SAR / head | < 50 % | Reconstruction matrix | 96 | Water suppression | no |
| Whole body / level | < 0.1 W/kg / | SENSE | yes | MTC | no |
| | normal | P reduction (AP) | 2.8 | Research prepulse | no |
| B1 rms | 0.93 uT | P os factor | 1 | Diffusion mode | no |
| PNS / level // VUIIS : | 59 % / normal | k-t BLAST | no | Elastography mode | no |
| dortch : | | Stacks | 1 | SAR mode | low |
| Sound Pressure Level | 28.01657 | type | parallel | B1 mode | default |
| (dB) | | slices | 46 | SAR Patient data | auto |
| MOTIO | N | slice gap | user defined | PNS mode | low |
| Cardiac synchronization | no | gap (mm) | 0 | 7 | |
| Heart rate > 250 bpm | no | slice orientation | transverse | Gradient mode | full control |
| Respiratory | no | l} | AP | max strength | 33 |
| compensation | | fold-over direction | | (mT/m) | 400 |
| Navigator respiratory | no | fat shift direction | Р | max slew rate (T/m/s) | 130 |
| comp | | Minimum number of | 1 | (1/111/5) | |
| Flow compensation | no | packages | 1.0.1 | | |
| Temporal slice spacing | equidistant | Slice scan order | default | _ | |
| fMRI echo stabilisation | no | PlanAlign | no | | |
| NSA | 1 | REST slabs | 0 | | |
| | | Catheter tracking | no | | |
| DYN/AN | | Interactive positioning | no | | |
| Angio / Contrast enh. | no | Allow table movement | no |] | |
| Quantitative flow | no | OFFC/AN | iG | 1 | |
| Manual start | yes | Stacks | 1 | 1 | |
| +Abuse dynamic loop | no | Stack Offc. AP | -5.532147 | - | |
| Dynamic study | individual | (P=+mm) | -3.332147 | | |
| dyn scans | 250 | RL (L=+mm) | -3.434925 | 1 | |
| recon multiplier | 1 | FH (H=+mm) | 35.45897 | 1 | |
| dyn scan times | shortest | Ang. AP (deg) | 0 | - | |
| FOV time mode | | | | | |
| | default | | - | 1 | |
| | default | RL (deg) | -5.457352 | | |
| dummy scans | 5 | | - | | |
| dummy scans immediate | | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction | 5 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan | 5 no no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device | 5 no no yes | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. | no no yes | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) | 5 no no yes 1 1119 | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization | 5 no no yes 1 1119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion | 5 no no yes 1 1119 | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. | 5 no no yes 1 1119 no no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole | 5 no no yes 1 1119 no no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling | 5 no no yes 1 1119 no no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRe | 5 no no yes 1 1119 no no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling | 5 no no yes 1 1119 no no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRe | 5 no no yes 1 119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO | 5 no no yes 1 119 no no no no no no no no full | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRt Preparation phases Interactive FO SENSE ref. scan | 5 no no yes 1 1119 no no no OC full no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey | 5 no no yes 1 1119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRt Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map | 5 no so yes 1 119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SEMSE ref. scan SmartPlan survey B0 field map B1 field map | 5 no no yes 1 1119 no no no OC full no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR | 5 no no yes 1 1119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map MIP/MPR Images | 5 no no yes 1 1119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PR: Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image | 5 no solution no s | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images | 5 no no yes 1 119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | 5 no no yes 1 1119 no no no OC full no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images | 5 no no yes 1 1119 no no no OC full no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | 5 no no yes 1 1119 no no no OC full no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey BO field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction | 5 no no yes 1 1119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PR: Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode | 5 no no yes 1 119 no no no no no no M, no, no, no, no Grey matter no soft | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PR! Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory | 5 no soft real time | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey BO field map B1 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data | 5 no no yes 1 1119 no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data Hardcopy protocol | 5 no no yes 1 119 no no no no no no M, no, no, no, no Grey matter no soft real time no no no no | RL (deg) | -5.457352 | | |
| dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey BO field map B1 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data | 5 no no yes 1 1119 no | RL (deg) | -5.457352 | | |


| | 41021 CEST fMRI | (13) 52:37.4 FMRI_ | TRAIL MAKING 04 | 1·14 በ | | |
|---|--|-----------------------------------|----------------------------|----------------------|--------------|--|
| INFO PAGE | | GEOMETRY | | CONTRAST | | |
| Total scan duration | 04:14.0 | Multi-transmit | no | Scan type | Imaging | |
| | 100 | Nucleus | H1 | Scan mode | MS | |
| Rel. signal level (%) Act. TR/TE (ms) | 2000 / 25 | Coil selection 1 | RX-Intf-1 | | FFE | |
| | | | | technique | | |
| Dyn. scan time | 00:02.000 | Xmit Coil selection | MTX-Volume- T/R | + ZOOM | no | |
| Time to k0 | 00:15.0 | User def elem sel | no | Contrast enhancement | no | |
| ACQ matrix M x P | 96 x 95 | element selection | All | Acquisition mode | cartesian | |
| ACQ voxel MPS (mm) | 2.50 / 2.50 / 2.50 | | | Fast Imaging mode | EPI | |
| DEC MDC (man) | | connection | conn-A | shot mode | single-shot | |
| REC voxel MPS (mm) | 2.50 / 2.50 / 2.50 | Coil selection 2 | RX-Intf-2 | Echoes | 1 | |
| Scan percentage (%) | 100 | element selection | All | partial echo | no | |
| | 1 | Dual coil | yes | shifted echo | no | |
| Packages | | CLEAR | yes | TE | user defined | |
| Min. slice gap (mm) | 0 | body tuned | yes | (ms) | 25 | |
| EPI factor | 37 | FOV RL (mm) | 240 | Flip angle (deg) | 63 | |
| Act. WFS (pix) / BW | 21.825 / 46.4 | AP (mm) | 240 | TR | user defined | |
| (Hz) | | FH (mm) | 115 | (ms) | 2000 | |
| BW in EPI freq. dir. (Hz) | 2878.9 | Voxel size RL (mm) | 2.5 | Halfscan | no | |
| Min. WFS (pix) / Max. | 21.770 / 46.5 | AP (mm) | 2.5 | Water-fat shift | minimum | |
| BW (Hz) | | Slice thickness (mm) | 2.5 | Shim | auto | |
| Min. TR/TE (ms) | 1999 / 12 | Recon voxel size (mm) | 2.5 | mDIXON | | |
| RF avg power computed | 1.371604 | Fold-over suppression | no | l | no | |
| (W) | F0 0: | Reconstruction matrix | 96 | Fat suppression | no | |
| SAR / head | < 50 % | SENSE SENSE | | Water suppression | no | |
| Whole body / level | < 0.1 W/kg / | l————————— | yes | MTC | no | |
| | normal | P reduction (AP) | 2.8 | Research prepulse | no | |
| B1 rms | 0.93 uT | P os factor | 1 | Diffusion mode | no | |
| PNS / level // VUIIS : | 59 % / normal | k-t BLAST | no | Elastography mode | no | |
| dortch : | | Stacks | 1 | SAR mode | low | |
| Sound Pressure Level | 28.01657 | type | parallel | B1 mode | default | |
| (dB) | | slices | 46 | SAR Patient data | auto | |
| MOTION | | slice gap | user defined | PNS mode | low | |
| Cardiac synchronization | no | gap (mm) | 0 | Gradient mode | full control | |
| Heart rate > 250 bpm | no | slice orientation | transverse | max strength | 33 | |
| Respiratory | no | fold-over direction | AP | (mT/m) | " | |
| compensation | | fat shift direction | Р | max slew rate | 130 | |
| Navigator respiratory | no | Minimum number of | 1 | (T/m/s) | | |
| comp | | packages | ' | <u> </u> | | |
| Flow compensation | no | Slice scan order | default | ł | | |
| Temporal slice spacing | equidistant | PlanAlign | no | • | | |
| fMRI echo stabilisation | no | REST slabs | 0 | · | | |
| NSA | 1 | | - | | | |
| DYN/AN | G | Catheter tracking | no | | | |
| Angio / Contrast enh. | no | Interactive positioning | no | | | |
| Quantitative flow | no | Allow table movement | no | | | |
| Manual start | | OFFC/AN | IG | | | |
| +Abuse dynamic loop | yes no | Stacks | 1 | | | |
| | individual | Stack Offc. AP | -5.532147 | | | |
| | | (P=+mm) | | l | | |
| Dynamic study | | | | ł | | |
| dyn scans | 120 | RL (L=+mm) | -3.434925 | | | |
| dyn scans recon multiplier | 120 | RL (L=+mm) FH (H=+mm) | -3.434925 35.45897 | | | |
| dyn scans recon multiplier dyn scan times | 120 1 shortest | | | | | |
| dyn scans recon multiplier dyn scan times FOV time mode | 120 1 shortest default | FH (H=+mm) | 35.45897 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans | 120 1 shortest default 5 | FH (H=+mm) Ang. AP (deg) | 35.45897 0 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate | 120 1 shortest default | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction | 120 1 shortest default 5 | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan | 120 1 shortest default 5 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction | 120 1 shortest default 5 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan | 120 1 shortest default 5 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device | 120 1 shortest default 5 no no yes | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. | 120 1 shortest default 5 no ves 1 | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion | 120 1 shortest default 5 no no yes 1 1119 | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization | 120 1 shortest default 5 no no yes 1 1119 | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole | 120 1 shortest default 5 no no yes 1 1119 | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. | 120 1 shortest default 5 5 no no yes 1 1119 no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole | 120 1 shortest default 5 no no yes 1 1119 no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling | 120 1 shortest default 5 no no yes 1 1119 no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC | 120 1 shortest default 5 5 no no yes 1 119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO | 120 1 shortest default 5 5 no no yes 1 1119 no no no OC full no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan | 120 1 shortest default 5 no no yes 1 1119 no no no fo fo full no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey | 120 1 shortest default 5 5 no no no yes 1 1 119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map | 120 1 shortest default 5 5 no no yes 1 1119 no no no DC full no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map | 120 1 shortest default 5 5 no no no yes 1 1119 no no no DC full no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR | 120 1 shortest default 5 5 no no no yes 1 1119 no no no DC full no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images | 120 1 shortest default 5 no no yes 1 1119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image | 120 1 shortest default 5 5 no no yes 1 1119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images | 120 1 shortest default 5 no no yes 1 1119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image | 120 1 shortest default 5 5 no no yes 1 1119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images | 120 1 shortest default 5 5 no no no yes 1 1119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | 120 1 shortest default 5 no no yes 1 119 no no no OC full no no no no M, no, no, no, no Grey matter | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast | 120 1 shortest default 5 5 no no yes 1 1119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode | 120 1 shortest default 5 5 no no yes 1 1119 no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SEINSE ref. scan SmartPlan survey B0 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory | 120 1 shortest default 5 no no no yes 1 119 no no no no o OC full no no no no M, no, no, no M no, no, no, no Grey matter no soft real time no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRO Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data | 120 1 shortest default 5 no no no yes 1 119 no no no OC full no no no M, no, no, no M, no, no, no Grey matter no soft real time no no no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRC Preparation phases Interactive FO SENSE Fef. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data Hardcopy protocol | 120 1 shortest default 5 no no yes 1 119 no no no no no M, no, no, no, no M no, no, no, no Grey matter no soft real time no no no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |
| dyn scans recon multiplier dyn scan times FOV time mode dummy scans immediate subtraction fast next scan synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRO Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data | 120 1 shortest default 5 no no no yes 1 119 no no no OC full no no no M, no, no, no M, no, no, no Grey matter no soft real time no no no no | FH (H=+mm) Ang. AP (deg) RL (deg) | 35.45897 0 -5.457352 | | | |

| | 41021 CEST fMRI | (13) 52:37.4 FMRI_ | nback 08:30.0 | | | |
|--|---|-------------------------|---------------|----------------------|--------------|--|
| INFO PAGE | | GEOMETRY | | CONTRAST | | |
| Total scan duration 08:30.0 | | Multi-transmit no | | Scan type | Imaging | |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | MS | |
| Act. TR/TE (ms) | 2500 / 25 | Coil selection 1 | RX-Intf-1 | technique | FFE | |
| Dyn. scan time | 00:02.500 | Xmit Coil selection | MTX-Volume- | + ZOOM | no | |
| Time to k0 | 00:18.7 | | T/R | Contrast enhancement | no | |
| ACQ matrix M x P | 96 x 95 | User def elem sel | no | Acquisition mode | cartesian | |
| ACQ voxel MPS (mm) | 2.50 / 2.50 / | element selection | All | Fast Imaging mode | EPI | |
| ` ′ | 2.50 | connection | conn-A | shot mode | single-shot | |
| REC voxel MPS (mm) | 2.50 / 2.50 / | Coil selection 2 | RX-Intf-2 | Echoes | 1 | |
| | 2.50 | element selection | All | partial echo | no | |
| Scan percentage (%) | 100 | Dual coil | yes | shifted echo | no | |
| Packages | 1 | CLEAR | yes | TE | user defined | |
| Min. slice gap (mm) | 0 | body tuned | yes | (ms) | 25 | |
| EPI factor | 37 | FOV RL (mm) | 240 | Flip angle (deg) | 63 | |
| Act. WFS (pix) / BW | 21.825 / 46.4 | AP (mm) | 240 | TR | user defined | |
| (Hz) | | FH (mm) | 115 | (ms) | 2500 | |
| BW in EPI freq. dir. (Hz) | 2878.9 | Voxel size RL (mm) | 2.5 | Halfscan | no | |
| Min. WFS (pix) / Max. | 21.770 / 46.5 | AP (mm) | 2.5 | Water-fat shift | minimum | |
| BW (Hz) | | Slice thickness (mm) | 2.5 | Shim | auto | |
| Min. TR/TE (ms) | 1999 / 12 | Recon voxel size (mm) | 2.5 | mDIXON | | |
| RF avg power computed | 1.097284 | Fold-over suppression | no | 11 | no | |
| (W) | 40.04 | Reconstruction matrix | 96 | Fat suppression | no | |
| SAR / head | < 40 % | SENSE SENSE | | Water suppression | no | |
| Whole body / level | < 0.1 W/kg / | P reduction (AP) | yes 2.8 | MTC | no | |
| D1 | normal | · · · · · · | | Research prepulse | no | |
| B1 rms | 0.84 uT | P os factor | 1 | Diffusion mode | no | |
| PNS / level // VUIIS : | 59 % / normal | k-t BLAST | no | Elastography mode | no | |
| dortch : | 27 12000 | Stacks | 1 | SAR mode | low | |
| Sound Pressure Level (dB) | 27.13099 | type | parallel | B1 mode | default | |
| MOTIO | u M | slices | 46 | SAR Patient data | auto | |
| | | slice gap | user defined | PNS mode | low | |
| Cardiac synchronization | no | gap (mm) | 0 | Gradient mode | full control | |
| Heart rate > 250 bpm | no | slice orientation | transverse | max strength | 33 | |
| Respiratory | no | fold-over direction | AP | (mT/m) | | |
| compensation Navigator respiratory | no | fat shift direction | Р | max slew rate | 130 | |
| comp | 110 | Minimum number of | 1 | (T/m/s) | | |
| Flow compensation | no | packages | | | | |
| Temporal slice spacing | equidistant | Slice scan order | default | | | |
| | · · | PlanAlign | no | | | |
| fMRI echo stabilisation | no | REST slabs | 0 | | | |
| NSA | 1 | Catheter tracking | no |] | | |
| DYN/AN | | Interactive positioning | no |] | | |
| Angio / Contrast enh. | no | Allow table movement | no | 1 | | |
| Quantitative flow | no | OFFC/AN | iG | 1 | | |
| Manual start | yes | Stacks | 1 | | | |
| +Abuse dynamic loop | no | Stack Offc. AP | -5.532147 | | | |
| Dynamic study | individual | (P=+mm) | 0.002117 | | | |
| dyn scans | 197 | RL (L=+mm) | -3.434925 | İ | | |
| recon multiplier | 1 | FH (H=+mm) | 35.45897 | | | |
| dyn scan times | shortest | Ang. AP (deg) | 0 | | | |
| FOV time mode | default | RL (deg) | -5.457352 | 1 | | |
| dummy scans | 5 | FH (deg) | 0 | 1 | | |
| immediate | no | (409) | 1- | 1 | | |
| subtraction | 1 | | | ~ | | |
| | | | | • | | |
| fast next scan | no | | | | | |
| fast next scan synch. ext. device | no yes | | | | | |
| synch. ext. device start at dyn. | | | | | | |
| synch. ext. device start at dyn. interval (dyn) | yes | | | | | |
| synch. ext. device start at dyn. | yes 1 | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion | yes 1 119 | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. | yes 1 1119 no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole | yes 1 1 1119 no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling | yes 1 1119 no no no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling | yes 1 1119 no no no no OC | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling | yes 1 1119 no no no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling | yes 1 1119 no no no no OC | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases | yes 1 119 no no no coc full | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 | yes 1 119 no no no o o o full no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan | yes 1 119 no no no co full no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey | yes 1 119 no no no CC full no no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRt Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map | yes 1 119 no no no coc full no no no no no no no no no no no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SEMSE ref. scan SmartPlan survey B0 field map B1 field map | yes 1 119 no no no no no no no no no no no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRi Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map MIP/MPR Images | yes 1 119 no no no no no no no no no no no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PR: Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image | yes 1 119 no no no no C full no no no no M, no, no, no, no M | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Impages Autoview image Calculated images | yes 1 119 no no no o C full no no no no M, no, no, no, no, no no, no, no, no, no, no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue | yes 1 119 no no no OC full no no no M, no, no, no, no Grey matter | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction | yes 1 119 no no no no OC full no no no no M, no, no, no M, no, no, no Grey matter no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PR Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast | yes 1 119 no no no no C full no no no M, no, no, no, no Grey matter no soft | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode | yes 1 119 no no no no CC full no no no no no M, no, no, no, no Grey matter no soft real time | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PR! Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory | yes 1 119 no no no no OC full no no no no M, no, no, no M no, no, no, no Grey matter no soft real time no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data | yes 1 119 no no no oc full no no no no M, no, no, no, no Grey matter no soft real time no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive FO SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data Hardcopy protocol | yes 1 119 no no no no CC full no no no M, no, no, no, no Grey matter no soft real time no no no | | | | | |
| synch. ext. device start at dyn. interval (dyn) dyn stabilization prospect. motion corr. Keyhole Arterial Spin labeling POST/PRI Preparation phases Interactive F0 SENSE ref. scan SmartPlan survey B0 field map B1 field map MIP/MPR Images Autoview image Calculated images Reference tissue EPI 2D phase correction Preset window contrast Reconstruction mode reuse memory Save raw data | yes 1 119 no no no oc full no no no no M, no, no, no, no Grey matter no soft real time no no | | | | | |


| | | (13) 52:37.4 🖵 T1_3D | _TFE_quantGeo | 00:55.7 | |
|------------------------------------|----------------------|-------------------------|---------------|----------------------|--------------|
| INFO PAGE | | GEOMETRY | | CONTRA | ST |
| Total scan duration | 00:55.7 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | 3D |
| Act. TR/TE (ms) | 2.8 / 1.44 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| ACQ matrix M x P | 192 x 190 | Xmit Coil selection | MTX-Volume- | + ZOOM | no |
| ACQ voxel MPS (mm) | 1.25 / 1.26 / | | T/R | Contrast enhancement | T1 |
| | 2.50 | User def elem sel | no | Acquisition mode | cartesian |
| REC voxel MPS (mm) | 0.94 / 0.94 / | element selection | All | Fast Imaging mode | TFE |
| | 2.50 | connection | conn-A | 3D non-selective | no |
| Scan percentage (%) | 98.9899 | Coil selection 2 | RX-Intf-2 | shot mode | multishot |
| TFE shots | 13 | element selection | All | TFE factor | 256 |
| TFE dur. shot / acq (ms) | 1675.7 / 729.0 | Dual coil | yes | 3D free factor | no |
| Min. TI delay | 386.7124 | CLEAR | yes | startup echoes | default |
| Act. WFS (pix) / BW | 0.778 / 1302.1 | body tuned | yes | +TFE followup echoes | 0 |
| (Hz) | / / / | FOV RL (mm) | 240 | shot interval | user defined |
| Min. WFS (pix) / Max. | 0.774 / 1308.4 | AP (mm) | 240 | (ms) | 4500 |
| BW (Hz) | 0.0000400 | FH (mm) | 165 | profile order | linear |
| RF avg power computed (W) | 0.8993402 | Voxel size RL (mm) | 1.25 | turbo direction | radial |
| · · | . 22.0/ | AP (mm) | 1.254902 | CENTRA (spiral) | no |
| SAR / head | < 33 % | FH (mm) | 2.5 | | |
| Whole body / level | 0.0 W/kg / normal | Recon voxel size (mm) | 0.94 | Echoes | 1 |
| D1 rmc | 0.76 uT | Fold-over suppression | no | partial echo | no |
| B1 rms | | | | shifted echo | no |
| PNS / level // VUIIS : dortch : | 59 % / normal | Slice oversampling | default | TE | shortest |
| Sound Pressure Level | 20 70210 | RF select. FOS | no | Flip angle (deg) | 7 |
| (dB) | 29.78219 | Reconstruction matrix | 256 | TR | shortest |
| MOTIO | \ | SENSE | yes | Halfscan | no |
| | - | P reduction (AP) | 2 | Water-fat shift | minimum |
| Cardiac synchronization | no | P os factor | 1 | Shim | auto |
| Heart rate > 250 bpm | no | S reduction (FH) | 2 | mDIXON | no |
| Respiratory | no | k-t BLAST | no | Fat suppression | no |
| compensation | | Overcontiguous slices | no | Water suppression | no |
| Navigator respiratory | no | Stacks | 1 | TFE prepulse | invert |
| comp | | slices | 66 | slice selection | no |
| Flow compensation | no | slice orientation | transverse | delay | user defined |
| fMRI echo stabilisation | no | fold-over direction | AP | (ms) | 1300 |
| Motion smoothing | no | fat shift direction | L | PSIR | no |
| NSA | 1 | Chunks | 1 | +inv pulse type | +B1 opt (low |
| DYN/AN | | PlanAlign | no | +iiiv puise type | BW) |
| Angio / Contrast enh. | no | REST slabs | 0 | мтс | no |
| Quantitative flow | no | Catheter tracking | no | T2prep | no |
| CENTRA | no | Interactive positioning | no | Research prepulse | no |
| Manual start | no | Allow table movement | no | Diffusion mode | no |
| + Abuse dynamic loop | no | | | | - |
| Dynamic study | no | OFFC/AN | | Elastography mode | no high |
| Arterial Spin labeling | no | Stacks | 1 | SAR mode | high |
| POST/PRO | oc | Stack Offc. AP | -4.247866 | B1 mode | default |
| Preparation phases | auto | (P=+mm) | 2.424025 | SAR Patient data | auto |
| Interactive F0 | no | RL (L=+mm) | -3.434925 | PNS mode | low |
| SENSE ref. scan | no | FH (H=+mm) | 35.88707 | Gradient mode | full control |
| SmartPlan survey | no | Ang. AP (deg) | 0 | max strength | 33 |
| B0 field map | no | RL (deg) | -5.457352 | (mT/m) | |
| B1 field map | no | FH (deg) | 0 | max slew rate | 166 |
| MIP/MPR | no | 1 | | (T/m/s) | <u> </u> |
| | | 1 | | | |
| Images Autovious image | M, no, no, no | 1 | | | |
| Autoview image | M | | | | |
| Calculated images | no, no, no, no | 1 | | | |
| Reference tissue | Grey matter | | | | |
| Preset window contrast | soft | | | | |
| B 1 1 1 1 1 | immediate | I | | | |
| Reconstruction mode | | 1 | | | |
| Save raw data | no | | | | |
| Save raw data Hardcopy protocol | no no | | | | |
| Save raw data | no | | | | |
| Save raw data Hardcopy protocol | no no | | | | |

| 🗀 Hospital (2) 📙 201 | 41021 CEST fMRI | (13) 52:37.4 T2star | _multiEcho 01:02 | 2.7 | |
|---|-----------------------|-------------------------|------------------|----------------------|--------------|
| INFO PAGE | | GEOMETRY | | CONTRAST | |
| Total scan duration | 01:02.7 | Multi-transmit | no | Scan type | Imaging |
| Rel. signal level (%) | 100 | Nucleus | H1 | Scan mode | 3D |
| Act. TR/TE1/delta TE | 34 / 3.3 / 3.2 | Coil selection 1 | RX-Intf-1 | technique | FFE |
| (ms) | | Xmit Coil selection | MTX-Volume- | loop order | zy_order |
| ACQ matrix M x P | 240 x 240 | | T/R | + ZOOM | no |
| ACQ voxel MPS (mm) | 1.00 / 1.00 / | User def elem sel | no | Contrast enhancement | T1 |
| 550 1100 () | 5.00 | element selection | All | Acquisition mode | cartesian |
| REC voxel MPS (mm) | 0.94 / 0.94 / 5.00 | connection | conn-A | Fast Imaging mode | none |
| Scan percentage (%) | 100 | Coil selection 2 | RX-Intf-2 | 3D non-selective no | |
| Act. WFS (pix) / BW | 1.418 / 714.8 | element selection | All | Echoes | 10 |
| (Hz) | 1.410 / /14.0 | Dual coil | yes | partial echo | no |
| Min. WFS (pix) / Max. | 1.125 / 900.9 | CLEAR | yes | shifted echo | no |
| BW (Hz) | | body tuned | yes | TE first | shortest |
| RF avg power computed | 0.1353424 | FOV AP (mm) | 240 | echospacing | shortest |
| (W) . | | RL (mm) | 240 | flyback | yes |
| SAR / head | < 5 % | FH (mm) | 60 | Flip angle (deg) | 8 |
| Whole body / level | 0.0 W/kg / | Voxel size AP (mm) | 1 | TR | shortest |
| | normal | RL (mm) | 1 | Halfscan | no |
| B1 rms | 0.29 uT | FH (mm) | 5 | Water-fat shift | user defined |
| PNS / level // VUIIS : | 60 % / normal | Recon voxel size (mm) | 0.9375 | (pixels) | 1.4 |
| dortch : | | Fold-over suppression | no | Shim | PB-volume |
| Sound Pressure Level | 29.7646 | Slice oversampling | default | ShimAlign | no |
| (dB) | | RF select. FOS | no | mDIXON | no |
| MOTIO | 1 | Reconstruction matrix | 256 | Fat suppression | no |
| Cardiac synchronization | no | SENSE | yes | Water suppression | no |
| Heart rate > 250 bpm | no | P reduction (RL) | 2 | MTC | no |
| Respiratory compensation | no | P os factor | 1 | Research prepulse | no |
| | no | S reduction (FH) | 1 | Diffusion mode | no |
| Navigator respiratory comp | no | k-t BLAST | no | Elastography mode | no |
| Flow compensation | yes | Overcontiguous slices | no | SAR mode | low |
| fMRI echo stabilisation | no | Stacks | 1 | B1 mode | default |
| NSA | 1 | slices | 12 | SAR Patient data | auto |
| DYN/AN | | slice orientation | transverse | PNS mode | low |
| Angio / Contrast enh. | inflow | fold-over direction | RL | Gradient mode | maximum |
| Quantitative flow | no | fat shift direction | Р | SofTone mode | no |
| Tone pulse | no | Chunks | 1 | Sol Folic Hiode | 110 |
| Manual start | no | PlanAlign | no | | |
| +Abuse dynamic loop | no | REST slabs | 0 | | |
| Dynamic study | no | Catheter tracking | no | | |
| Arterial Spin labeling | no | Interactive positioning | no | | |
| POST/PR | | Allow table movement | no | | |
| Preparation phases | auto | OFFC/ANG | | | |
| _ <u>' </u> | | Stacks | 1 | | |
| Interactive F0 SENSE ref. scan | no no | Stack Offc. AP | -5.532147 | | |
| SmartPlan survey | no | (P=+mm) | | | |
| B0 field map | | RL (L=+mm) | -3.434925 | | |
| B1 field map | no no | FH (H=+mm) | 35.45897 | | |
| MIP/MPR | | Ang. AP (deg) | 0 | | |
| | no M, R, I, no | RL (deg) | -5.457352 | | |
| Images Autoview image | M, R, I, no | FH (deg) | 0 | | |
| Calculated images | T2, no, no, no | Shim Size AP (mm) | 180.7495 | | |
| T2* clipvalue (ms) | 12, no, no, no | RL (mm) | 134.1432 | | |
| Reference tissue | Grey matter | FH (mm) | 114.2866 | | |
| | soft | Offc. AP (P=+mm) | -1.87074 | | |
| Preset window contrast | | RL (L=+mm) | -3.660397 | | |
| Reconstruction mode | immediate | FH (H=+mm) | 33.86363 | | |
| Save raw data | no | Ang. AP (deg) | 2.26406 | | |
| Hardcopy protocol | no | RL (deg) | -6.580403 | | |
| Ringing filtering | default | FH (deg) | 4.637685 | | |
| Geometry correction | default | , <i>y</i> | | i . | |
| Elliptical k-space shutter | | ŧ | | | |

Exam Card

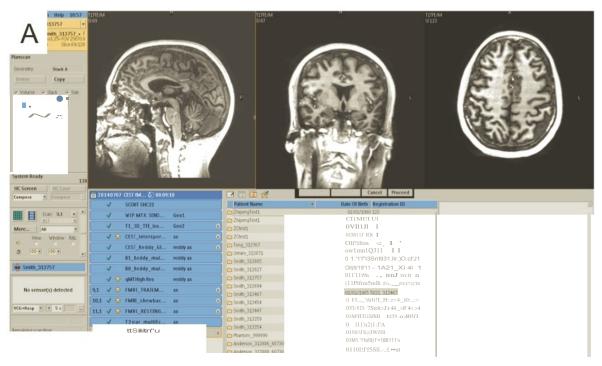

Viewing Environment

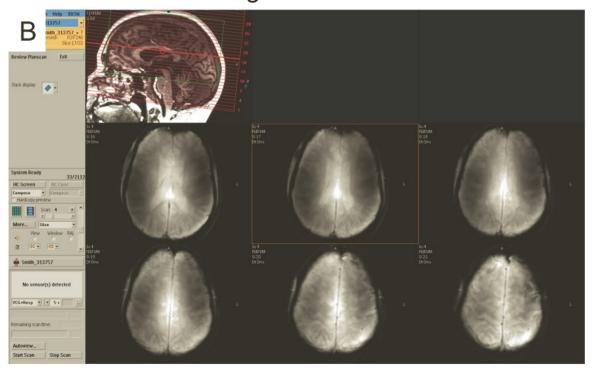



Contents lists available at SciVerse ScienceDirect

Neurolmage

journal homepage: www.elsevier.com/locate/ynimg



Quantitative magnetization transfer imaging of human brain at 7 T

Richard D. Dortch ^{a,b,*}, Jay Moore ^{a,b}, Ke Li ^{a,b}, Marcin Jankiewicz ^{a,b}, Daniel F. Gochberg ^{a,b,c}, Jane A. Hirtle ^b, John C. Gore ^{a,b,c,d,e}, Seth A. Smith ^{a,b,d}

- ^a Department of Radiology and Radiological Sciences, Vanderbilt University, Nashville, TN, USA
- ^b Vanderbilt University Institute of Imaging Science, Vanderbilt University, Nashville, TN, USA
- ^c Department of Physics and Astronomy, Vanderbilt University, Nashville, TN, USA
- ^d Department of Biomedical Engineering, Vanderbilt University, Nashville, TN, USA
- ^e Department of Molecular Physiology and Biophysics, Vanderbilt University, Nashville, TN, USA

article info

Article history:
Accepted 18 August 2012
Available online 23 August 2012

Keywords: Magnetization transfer 7 T White matter Myelin Brain Multiple sclerosis

abstract

Quantitative magnetization transfer (qMT) imaging yields indices describing the interactions between free water protons and immobile macromolecular protons. These indices include the macromolecular to free pool size ratio (PSR), which has been shown to be correlated with myelin content in white matter. Because of the long scan times required for whole-brain imaging (≈ 20-30 min), qMT studies of the human brain have not found widespread application. Herein, we investigated whether the increased signal-to-noise ratio available at 7.0 T could be used to reduce qMT scan times. More specifically, we developed a selective inversion recovery (SIR) qMT imaging protocol with a i) novel transmit radiofrequency (B_1^t) and static field (B_0) insensitive inversion pulse, ii) turbo field-echo readout, and iii) reduced TR. In vivo qMT data were obtained in the brains of healthy volunteers at 7.0 T using the resulting protocol (scan time ≈ 40 s/slice, resolution = 2 × 2 × 3 mm³). Reliability was also assessed in repeated acquisitions. The results of this study demonstrate that SIR qMT imaging can be reliably performed within the radiofrequency power restrictions present at 7.0 T, even in the presence of large B+ and B inhomogeneities. Consistent with qMT studies at lower field strengths, the observed PSR values were higher in white matter (mean \pm SD = 17.6 \pm 1.3%) relative to gray matter (10.3 \pm 1.6%) at 7.0 T. In addition, regional variations in PSR were observed in white matter. Together, these results suggest that qMT measurements are feasible at 7.0 T and may eventually allow for the high-resolution assessment of changes in composition throughout the normal and diseased human brain in vivo

© 2012 Elsevier Inc. All rights reserved.

Introduction

In addition to the free water protons typically observed in magnetic resonance imaging (MRI), there are protons residing on immobile macromolecules in tissue (Wolff and Balaban, 1989). Typical imaging sequences do not directly detect this pool of protons because they exhibit very short transverse relaxation times ($\approx 10~\text{lJs}$) and, therefore, lose coherence before their signal can be captured. This macromolecule proton pool can, however, be indirectly detected by exploiting its interactions with the free water pool via chemical exchange and/or dipolar mechanisms [referred to together as the magnetization transfer (MT) effect]. Previous phantom studies (Koenig, 1991; Kucharczyk et al., 1994) have shown that the bulk of the MT effect in white matter (WM) arises from myelin-associated lipids, which suggests that MT contrast may be

E-mail address: richard.dortch@vanderbilt.edu (R.D. Dortch).

a more specific marker for myelin pathology than conventional imaging methods. As a result, there is considerable interest in exploiting MT contrast to assay changes in myelination associated with a number of diseases [e.g., multiple sclerosis (Catalaa et al., 2000; Filippi and Rocca, 2004; Gass et al., 1994; Kalkers et al., 2001) and neuropsychiatric diseases (Bruno et al., 2004; Kabani et al., 2002a, 2002b)].

MT contrast can be generated by applying an off-resonance radiofrequency (RF) prepulse to selectively saturate the spectrally broad macromolecular proton pool (Wolff and Balaban, 1989). This saturation then transfers to the free water proton pool via MT, resulting in a decrease in the observed free water signal. The magnitude of this effect can be characterized by a semi-quantitative metric known as the magnetization transfer ratio (Dousset et al., 1992): $\text{MTR} = 1 - S_{\text{sat}}/S_0, \text{ where } S_{\text{sat}} \text{ and } S_0 \text{ are the observed signal intensities with and without the application of an MT saturation prepulse, respectively. Although the MTR has been shown to correlate with myelin content (Odrobina et al., 2005; Schmierer et al., 2004), it is also sensitive to the choice of experimental parameters such as RF power (Berry et al., 1999) as well as non-MT-specific NMR parameters such as tissue relaxation times (Henkelman et al., 1993). As a result, quantitative MT (gMT) approaches$

Grant sponsors: NIH K01 EB009120 (SAS), NIH T32 EB001628 (JCG), NIH EB00461 (JCG), and Vanderbilt Bridge Funding (DFG).

^{*} Corresponding author at: Vanderbilt University Institute of Imaging Science, AA-1101 Medical Center North, 1161 21st Avenue South, Nashville, TN 37232-2310, USA. Fax: +1 615 322 0734.

quantify distinct tissue characteristics (e.g., the size of the macromolecular pool, rate of MT exchange) rather than the combined effect of multiple tissue and/or acquisition parameters. As such, qMT measures are thought to yield more specific information on tissue composition than the MTR.

Pulsed saturation qMT imaging (Graham and Henkelman, 1997; Pike, 1996; Sled and Pike, 2000, 2001) has received considerable attention for application in humans in vivo because it allows for the rapid collection of gMT data within the hardware constraints of most clinical systems. This approach involves a steady-state, spoiled gradient-echo acquisition interleaved with an MT-preparation pulse. By collecting images over a range of MT pulse offset frequencies and/or powers and fitting the resulting data to a two-pool model of the MT effect, one can extract parameters such as the macromolecular to free pool size ratio (PSR) and the rate of MT exchange. Previous work has shown that the PSR is correlated with myelin content (Odrobina et al., 2005; Ou et al., 2009; Schmierer et al., 2007; Underhill et al., 2011). The relationship between the rate of MT exchange and underlying tissue composition is less clear; however, previous work has suggested that the rate of MT exchange may reflect changes within the myelin lipid structure (Smith et al., 2009).

Unfortunately, qMT imaging has not found widespread application in practice. This can be attributed in part to the long scan times (\approx 20–30 min for whole-brain imaging) required to collect images at multiple offset frequencies and/or powers. The number of total images required can be reduced by designing optimal sampling strategies (Cercignani and Alexander, 2006; Levesque et al., 2011) or by fixing certain model parameters in the fitting procedure (Underhill et al., 2009, 2011). Potentially more efficient strategies based upon steady-state free-precession (SSFP) sequences (Garcia et al., 2010; Gloor et al., 2008) may also be employed.

As an alternative, or perhaps in combination with these strategies, one could translate qMT imaging approaches to higher field strengths. The resulting increase in SNR could then be used to obtain more reliable estimates of MT parameters or traded to reduce scan times and/or increase resolution. To date, qMT studies in humans in vivo have been primarily limited to 1.5 and 3.0 T and we are aware of only one report (Mougin et al., 2010) of MT parameters in the human brain in vivo at 7.0 T. The translation of pulsed saturation approaches to 7.0 T faces two primary challenges: i) RF power limitations [e.g., specific absorption ratio (SAR) limitations] and ii) transmit RF (B^+) and static magnetic field (B_0) inhomogeneities. SSFP-based approaches may also be limited at high field by banding artifacts associated with B_0 inhomogeneities. In contrast, selective inversion recovery (SIR) qMT

imaging (Edzes and Samulski, 1977; Gochberg et al., 1997), which is based upon measuring the biexponential recovery of the free water pool in the presence of MT after an on-resonance inversion pulse, has been suggested (Dortch et al., 2011) to be less sensitive to these issues. Note that this approach is similar to the stimulated echo approach proposed by Ropele et al. (2003); therefore, both approaches may be well suited for qMT imaging at 7.0 T.

In this study, we have investigated the feasibility of using the SIR approach for high field qMT imaging of the human brain. More specifically, we have translated our previously published 3.0-T SIR protocol (Dortch et al., 2011) to 7.0 T with two significant modifications. First, we incorporated a novel B_1^+ - and !1 B_0 -insensitive composite inversion pulse to ensure a more uniform inversion of the free water pool over the whole brain. Second, we transitioned from a turbo-spin echo readout (TSE) to a turbo field-echo readout (TFE)—similar to an MP-RAGE sequence (Mugler and Brookeman, 1990)—as the former is susceptible to B_1^+ -related artifacts (due to imperfect refocusing) and is SAR-limited at high field. The TFE readout has the added benefit of covering k-space more efficiently than the TSE readout, which, in combination with some additional protocol optimization, allowed us to transition from a single-slice approach at 3.0 T to a whole-brain approach at 7.0 T (\approx 40 s/slice at 2.0 × 2.0 × 3.0 mm³ resolution). Using

this protocol, in vivo qMT data were obtained in the brains of 13 healthy volunteers at 7.0 T. To assess the reproducibility of the technique, six of the healthy volunteers were scanned twice. Additional numerical simulations were performed to determine the effect of TFE readout on our qMT parameter maps.

Theory

Consider free water (f) and macromolecular (m) proton pools between which MT can occur. Define unique equilibrium magnetizations (M_{0f} and M_{0m}), spin-lattice relaxation rates (R_{1f} and R_{1m}), and spin-spin relaxation rates (R_{2f} and R_{2m}) for each pool as well as an MT rate from the macromolecular to the free pool (k_{mf})—the rate in the other direction can be determined from $k_{fm} = k_{mf} M_{0m} / M_{0f}$. Assume MT of transverse magnetization to be negligible because of the short T_2 of the macromolecular pool. In this case, the transverse components of the macromolecular pool can be ignored. The time evolution of the remaining x, y, and z components of the magnetization vector $M = \left[M_{xf} \ M_{yf} \ M_{zf} \ M_{zm}\right]^T$ during a constant amplitude RF pulse can be expressed in matrix form as (Portnoy and Stanisz, 2007)

$$\frac{dM\delta t \flat}{dt} \frac{1}{4} AM\delta t \flat B; \qquad \qquad \delta 1 \flat$$

where

!1ω is the frequency offset from resonance for the RF pulse, $ω_1$ is the frequency of precession about the RF pulse, and φ is the phase of the RF pulse in the transverse plane. The standard Bloch equations implicitly assume a Lorentzian lineshape, which is invalid for the macromolecular proton pool. As a result, the Bloch equations for the macromolecular pool have been replaced in Eq. (2) by a single longitudinal component whose saturation is governed by the rate R_{RF} = T(α).

 2g (ω), where g is the lineshape function of the macromolecular pool. When applying off-resonance irradiation, a super-Lorentzian lineshape is typically used to model biological macromolecular protons (Morrison et al., 1995). Because the super-Lorentzian exhibits an on-resonance singularity, Gaussian (Gochberg and Gore, 2007) or super-Lorentzian functions extrapolated from a 1 kHz offset (Gloor et al., 2008) are typically used to model the macromolecular pool lineshape pool during on-resonance irradiation.

The general solution to this system of equations can be expressed as

Mỗt
$$\frac{1}{4}$$
 expổ A t $\frac{1}{4}$ exp δ At $\frac{1}{4}$

where M(0) is the initial condition of the system and I is an identity matrix. The same expression can be used to describe the system during free precession (i.e., when $\omega_1 = 0$). In this case, the solution can be further simplified by noting that the z-component is decoupled from the x- and y-components, resulting in the following expression for the longitudinal magnetization vector $M_z = [M_{zf} \ M_{zm}]^T$

where $M_0 = [M_{0f} M_{0m}]^T$ and A_z is the lower-right quadrant of A with $R_{RF} = 0$. Expanding the matrix exponentials in this expression yields

where $\lambda^{+\prime-}$ are the negative eigenvalues of A_z and U is a matrix whose columns are the corresponding eigenvectors. From Eq. (5), it can be seen that M_{zf} recovers as a biexponential function governed by the fast and slow rate constants λ^+ and λ^- , respectively, during free precession. As described below, one can obtain estimates of qMT parameters (e.g., PSR and $k_{mf})$ by measuring this biexponential recovery.

Methods

Pulse sequence

The SIR qMT sequence (Fig. 1) used herein is similar to the inversion recovery sequence used to measure T₁ with two modifications. First, short inversion times (≈ 10 ms or less) are sampled in order to capture the fast-recovering λ^+ component of the biexponential recovery. Second, a T2-selective inversion pulse is applied. This is achieved via a low power inversion pulse whose duration is much longer than the T_2 of the macromolecular pool ($T_{2m} \approx 10 \text{ Us}$) and much shorter than the T_2 of the free water pool ($T_{2f} \approx 10-100$ ms). Ideally, this pulse inverts M_{zf} with minimal saturation of M_{zm} . In other words, this pulse maximizes the difference between the pools and, in turn, the sensitivity of the signal to MT. This is followed by a variable duration inversion recovery period to sample the transient biexponential recovery of Mzf and a center-out TFE readout (SIR-TFE) to efficiently sample k-space. For inversion recovery acquisitions, a predelay time $\,t_{\text{d}} \approx \, 5/\lambda^{-}$ is commonly employed to ensure full recovery of Mzf. However, if one can assume that the longitudinal magnetization of both pools is approximately zero at the end of the readout, the effect of a shorter predelay period can be accounted for in the signal model, allowing one to reduce t_d (and scan times) without biasing the estimated parameters. This assumption has been previously shown to hold true for a TSE readout (Gochberg and Gore, 2007); however, this cannot be assumed for the TFE readout employed herein. As a result, we empirically designed a train of RF pulses [number of pulses = 32, α = 135°, pulse spacing = 20 ms, pulse train duration $(t_{sat}) = 620 \text{ ms}$ to saturate both pools following the TFE readout. To assess the effect of this pulse train on the longitudinal magnetization of both pools, numerical simulations were performed via Eq. (3) and the following parameters: $R_{1m} = R_{1f} = 0.8 \text{ s}^{-1}$, $T_{2m} = 10 \text{ !Js}$ (Gaussian

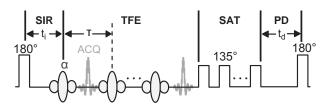


Fig. 1. SIR-TFE pulse sequence diagram. The sequence employs i) a composite inversion pulse (Fig. 2) designed to uniformly invert M_{zf} over a range of expected !1B0 and $B_1^{\ +}$ values with minimal macromolecular pool saturation, ii) a variable duration inversion recovery (SIR) period to sample the free pool recovery, iii) a TFE readout to efficiently cover k-space, iv) a pulse train to saturate (SAT) the free and macromolecular pools (allows t_d b5/ λ^-), and v) a predelay (PD) period to allow for partial M_z recovery. Legend: $t_1=$ inversion time, $t_d=$ predelay, $\tau=$ TFE pulse-to-pulse interval, ACQ = acquisition.

lineshape), $T_{2f}\!=\!60$ ms, $k_{mf}\!=\!15~s^{-1}$, and PSR = 15%. The results from these simulations indicate that this pulse train saturates both pools $[M_{zf}(t_{sat})/M_{0f} \leq 0.01$ and $M_{zm}(t_{sat})/M_{0m} \leq 0.06]$ over the range of expected B+ values (B+ /B+ = 0.3-1.0) in the human brain

1,actual 1,nominal

at 7.0 T [n.b., the manufacturer-provided power optimization tended to yield a mean $B_{1,actual}^+/B_{1,nominal}^+$ b 1.0 (Moore et al., 2010)].

Plugging the initial condition of $M_z(t_d\!=\!0)\!=\!0$ into Eq. (4), signal equations can be generated for the predelay period of the SIR-TFE sequence. The ending values for this period can then be used as the initial condition for the inversion recovery period, taking account for the effect of the inversion pulse

$$M_{z}\delta t^{b} p \% SM_{z}\delta t^{-} p;$$
 $\delta 6 p$

where S is a diagonal matrix with elements that account for the inversion of the free pool ($S_{\rm f}\!=\!-$ 1 denotes complete inversion) and the saturation of the macromolecular pool ($S_{\rm m}\!=\!1$ denotes no saturation) and $t^{+\prime-}$ is the time immediately before/after the pulse. This yields the final expression for the evolution of M_z during the SIR period of the sequence

$$M_z \delta t_i; t_d \triangleright \frac{1}{4} \delta \exp \delta A_z t_i \triangleright S \mathbb{I} - \exp \delta A_z t_i \triangleright \mathbb{I} - \exp \delta A_z t_i \triangleright \mathbb{I} - \exp \delta A_z t_i \triangleright S \mathbb{I}$$

In addition to these pulse sequence modifications, a novel 64-element composite inversion pulse was designed and employed herein to ensure a uniform inversion of $M_{\rm zf}$ over the range of $\mbox{\mbox{\it B}}^+$ and !1B $_{\rm 0}$ values previously measured in the human brain at 7.0 T (Moore et al., 2010). The optimization procedure (Moore et al., 2010) tended to produce high power pulses with suboptimal T $_{\rm 2}$ -selectivity. As a result, we included an additional RF power constraint into the procedure, which was weighted against the uniform inversion constraint. The resulting amplitudes and phases of the subpulses are shown in Fig. 2. To evaluate the pulse's performance, S $_{\rm f}$ and S $_{\rm m}$ were estimated from Eq. (6) by propagating Eq. (3) through each of the 64 subpulses [neglecting T $_{\rm 1}$ relaxation and MT during the pulse and

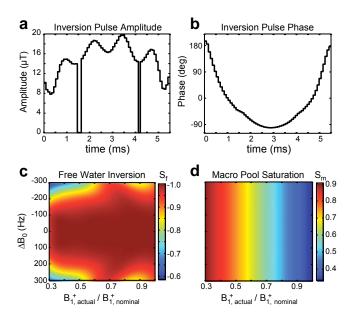


Fig. 2. Composite inversion pulse amplitudes (a), phases (b), predicted free water inversion efficiency $\S_t(c)$, and predicted macromolecular saturation fractions $\S_m(d)$. $\S_r = -1$ denotes complete inversion; $\S_m = 1$ denotes no saturation. The two zero-amplitude discontinuities in the RF pulse (a) are a consequence of the power constraint used in the minimization procedure. The RF phase (b) of the pulse at these discontinuities is arbitrary; therefore, the phase at these points was set based upon linear interpolation of the neighboring RF phases for display purposes.

assuming a Gaussian macromolecular pool lineshape with $T_{2m} = 10 \, \text{ls}$ (Gochberg and Gore, 2007)]. From this procedure, the pulse is predicted to yield a uniform inversion of M_{zf} over a wide range of B_1^+ and $!1B_0$ values without complete saturation of M_{zm} .

Numerical simulations

The TFE readout employed herein effectively blurs the image along the phase-encoding direction according to its readout point-spread function (PSF), which is a complex function of the sequence timings and the NMR parameters of the tissue (Constable and Gore, 1992). If the readout PSF is constant as a function of t_i , then its effect will be to simply blur the final MT parameter maps. If, however, the readout PSF changes as a function of t_i , each image will be blurred to a different degree, potentially biasing the final parameter maps.

To evaluate this effect, the SIR-TFE signal arising from a onedimensional (1D) test object was numerically simulated. As shown in Fig. 3, MT parameters were defined for test object regions representing white matter (WM), gray matter (GM), and cerebrospinal fluid (CSF). For each region and t_i, the signal evolution during each RF pulse and precession period of the TFE readout was simulated from Eq. (3) with the imaging parameters in the Data acquisition section—using $M_{z}(t_{i})$ from Eq. (7) as the initial condition and replacing each time-varying excitation pulse with a constant amplitude pulse of equivalent flip angle or root-mean-squared power (Ramani et al., 2002) for the free water or macromolecular pool, respectively. Complete spoiling of transverse magnetization was assumed prior to each RF pulse. The resulting $\ensuremath{\text{M}_{\text{zf}}}$ immediately after each RF pulse was taken to represent the signal as a function of echo number. The signal was then re-ordered to account for the k-space trajectory and SENSE acceleration used, and the resulting re-ordered signal was taken to represent a k-space filter. To apply the k-space filters to the 1D test object, each uniform object region was Fourier transformed into k-space, multiplied by its corresponding k-space filter, and inverse Fourier transformed back into image space. The resulting object regions were then summed to generate

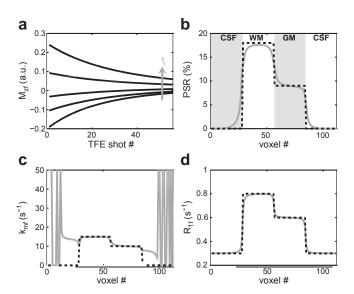


Fig. 3. Numerical simulations of the SIR-TFE readout (a) and resulting qMT parameter fits for the 1D test object defined in (b–d). (a) The $M_{\rm zf}$ for each region (WM shown here) and t_i was simulated and reordered into the corresponding k-space filter. The resulting filters were applied to the 1D test object as described in the text. (b–d) From the simulated fit parameters (solid gray lines), it can be seen that the TFE readout blurs the parameter maps with little or no bias (except for $k_{\rm mf}$ in CSF regions, which do not exhibit an MT effect).

the final blurred 1D object at each t_i . Finally, to assess the effect of the TFE readout on qMT parameter maps, the magnitude of the blurred test object signal at each voxel was fit to the M_{zf} component of Eq. (7) as described in the Data analysis section.

Subjects

MRI was performed on thirteen healthy volunteers (22–37 years old, 10 male, 3 female). To test reproducibility, six of the healthy volunteers were asked to undergo a second MRI scan at least two weeks after the first session. The study was approved by our local institutional review board, and signed consent was obtained prior to all examinations.

Data acquisition

Imaging was performed using a 7.0-T, Philips Achieva MR scanner (Philips Healthcare, Best, The Netherlands). A quadrature volume coil was used for excitation and a 32-channel head coil (Nova Medical, Wilmington, MA, USA) was used for signal reception. For qMT imaging, SIR-TFE data were collected in each subject using the general pulse sequence shown in Fig. 1.

An initial experiment was performed in one healthy volunteer to determine the effect of the post-TFE saturation train and predelay time $t_{\rm d}$ on the qMT parameter maps. For this initial experiment, SIR-TFE data were acquired in a single 5-mm axial slice with and without the post-TFE saturation train (see Fig. 1, shaded area labeled SAT) over a range of $t_{\rm d}$ values (0.125–10 s). Additional imaging parameters included: $t_{\rm i}$ logarithmically spaced between 6 ms and 2 s (15 values) and $t_{\rm i}$ = 10 s, TFE echoes per shot = 53, TFE pulse-to-pulse interval (τ)/TE/ α = 2.8 ms/1.4 ms/15°, SENSE factor = 2, field-of-view = 212 × 212 mm², resolution = 2.0 × 2.0 mm², and number of signal acquisitions averaged (NSA) = 2.

Based upon the results of this experiment along with previous numerical simulations (Gochberg and Gore, 2007), a $t_{\rm d}$ of 2.5 s was chosen to balance the scan time and SNR constraints for whole-brain SIR-TFE imaging. Whole-brain qMT data were acquired in 12 volunteers (six scanned twice) via a three-dimensional (3D) SIR-TFE sequence using the previously listed parameters except: $t_{\rm i}$ logarithmically spaced between 6 ms and 2 s (13 values) and $t_{\rm i}\!=\!8$ s, SENSE factor =4 (2 anterior–posterior, 2 superior–inferior), field-of-view =212×212×90 mm³, resolution =2.0×2.0×3.0 mm³, and NSA =1. This resulted in an acquisition time \approx 19 min for 30 slices.

Recall that the signal model [Eq. (7)] has terms (S_f and S_m) that account for the effect of the inversion pulse on the free and macromolecular pool magnetizations. S_f was included as a free parameter in the fit as described in the Data analysis section, while S_m was numerically estimated as described in the Pulse sequence section. Because S_m is sensitive to B_1^+ (see Fig. 2d), this numerical estimation required an independent measurement of B_1^+ . As a result, B_1^+ was estimated in same volume as the SIR-TFE data using the actual flip angle imaging (AFI) method (Yarnykh, 2007) with $TR_1/TR_2 = 125/25$ ms and a 60° slab-selective excitation pulse (asymmetric sinc pulse with Gaussian apodization).

Data analysis

All data analyses were performed in MATLAB (Mathworks, Natick, MA). Prior to data fitting, each SIR-TFE and AFI volume was co-registered to the SIR-TFE volume acquired at $t_i = 110$ ms (middle value) using a 3D rigid body registration based upon normalized mutual information (Viola and Wells, 1997). Following co-registration, automatic brain extraction was performed (Smith, 2002) and qMT parameter maps were calculated in each volunteer. The SIR-TFE signal model described in Eq. (7) has seven independent parameters: R_{1m} , R_{1f} , S_m , S_f , M_{0f} , $PSR = M_{0m}/M_{0f}$, and k_{mf} ($k_{fm} = k_{mf}PSR$). As is the

case with pulsed saturation methods, the signal dependence on R_{1m} for SIR data is weak (Li et al., 2010). Therefore, R_{1m} was set equal to R_{1f} for fitting purposes. The parameter S_m was numerically estimated for each voxel. This required an independent estimate of the actual flip angle in each voxel (α_{actual}), which was calculated from the AFI data using the following relationship (Yarnykh, 2007):

$$\alpha_{\text{actual}} \frac{1}{4} \cos^{-1} \frac{(n-1)}{n-r};$$
 õ8Þ

where n=TR₂/TR₁, r=S(TR₂)/S(TR₂), and S is the signal intensity. The resulting α_{actual} map was smoothed with a 10×10×9 mm³ moving-average filter to minimize the impact of imaging artifacts. Following this operation, the flip angle values were converted to B⁺_{1,actual} values for the composite inversion pulse (see Fig. 2a) and S_m was estimated using the procedure described in the Pulse sequence section (see Fig. 2d). The remaining five parameters (R_{1f}, S_f, M_{0f}, k_{mf}, and PSR) were estimated for each voxel by fitting SIR-TFE

data (14 $\,t_{\rm i}$ values) to the $\,M_{\rm zf}$ component of Eq. (7) in a least-squares sense using the procedure described in Dortch et al. (2011).

SIR-TFE data had a mean SNR per voxel of 180 ± 50 (range = 60-320) within the defined ROIs, where SNR is defined as M_{of} divided by the standard deviation (SD) of the residuals of the fit. Monte Carlo simulations, similar to those described by Li et al. (2010), were performed to predict the uncertainty of the fit parameters at these SNR levels. The t_i and t_d values listed above were used for these simulations. Additional simulation parameters included: $R_{1m}=R_{1f}=0.8~s^{-1}$, $k_{mf}=15~s^{-1}$, PSR = 15%, $s_m=0.7$, and $s_f=-0.95$. Over an

SNR range of 60–320, the SDs of the fit PSR, R_{1f} , and k_{mf} values were 0.4–2.2%, 0.01–0.03 s $^ ^1$, and 0.9–5.5 s $^ ^1$, respectively. This suggests

that PSR and R_{1f} can be robustly determined from the in vivo brain data collected herein. Consistent with previous studies (Li et al., 2010), the uncertainty in k_{mf} is expected to be much larger, especially in lower SNR regions.

Following this fitting procedure, qMT parameter maps were smoothed with a locally-adaptive Gaussian filter (kernel size = $10 \times 10 \times 9 \text{ mm}^3$, full width at half maximum = 1/2 kernel size) to remove outliers that tended to occur at tissue boundaries. To perform this operation, each filtered map was subtracted from the raw parameter map, and outliers were defined as voxels whose value was three standard deviations above the mean difference across all voxels. For these outliers, the value in the raw parameter map was replaced with the value in the filtered map. This process was iterated until the number of outliers was less than the expected value (0.3% of the total number of voxels).

Statistics

Mean qMT parameters (PSR, R_{1f} , and k_{mf}) were calculated within the following regions-of-interest (ROI): head of the caudate, putamen, thalamus, genu and splenium of the corpus callosum, internal capsule, corona radiata, occipital WM, and frontal WM. Statistical comparisons were performed on the mean ROI values to evaluate each parameter's i) variation across ROIs (i.e., regional differences), ii) variation and reproducibility across time, and iii) variation across volunteers. To compare parameters across WM regions, a non-parametric Wilcoxon rank-sum test was performed, with a pb 0.05 deeming a significant difference between ROI values. To evaluate the test-retest reproducibility of each parameter, a Bland-Altman (BA) analysis was performed. For the BA analysis, the mean difference and the limits of agreement (LOA = mean difference ± 1.96 * SD) were tabulated across scans for all ROIs. Additionally, a Wilcoxon signed-rank test was performed between the test and retest parameter values for each ROI, with a p > 0.05 indicating a non-significant difference between scans at each time point. To assess the test-retest variability of each parameter within each ROI, the coefficient of variation was calculated from: CV S= M

(2Þ 100, where S is the SD of the test–retest difference across subjects. M is the mean

s is the SD of the test-retest difference across subjects, M is the mean value across all test-retest scans and subjects, and the 2 term accounts for the propagation of uncertainty from the difference operation. The across-cohort variability of each parameter within each ROI was also assessed via: CV_{cohort} = S/M* 100, where S is the SD across the cohort and M is mean value across the cohort. All values are reported as the mean ± SD unless otherwise stated.

Results

The results of the numerical simulations designed to assess the effect of TFE readout on qMT parameter maps are shown in Fig. 3. In Fig. 3a, the evolution of $M_{\rm zf}$ during the TFE readout is shown for WM as a function of $t_{\rm i}$. Note that this evolution is related to the k-space filter of the readout. It can be seen that the shape (width and rate of decay) of the k-space filter changes as a function of $t_{\rm i}$, which manifests as a change in object blurring as a function of $t_{\rm i}$. The effect of this on the qMT parameter maps is shown in Figs. 3b–d. It can be seen that the resulting qMT parameter maps are smoothed in the phase-encoding direction with little bias in the fit parameters. It should be noted, however, that PSR values were slightly underestimated in the WM region of the 1D test object. Additional simulations indicated that this bias increased as the size of the WM region decreased

Fig. 4 displays PSR maps acquired with and without application of the post-TFE saturation train (see the SAT region in Fig. 1) as a function of $t_{\rm d}$. For scans with the saturation train, the fit PSR values at

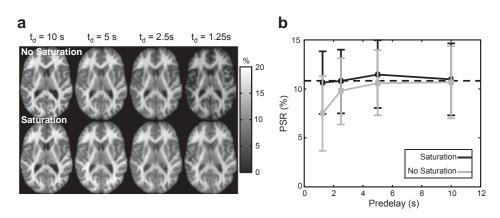


Fig. 4. (a) Maps of PSR as a function of t_d without (top) and with (bottom) a post-TFE saturation train and (b) corresponding mean (\pm SD) slice-wise PSR values. For scans with the saturation train, all PSR values were nearly identical to the values at full recovery (dashed line). Without the saturation train, PSR values were increasingly underestimated with decreasing t_d .

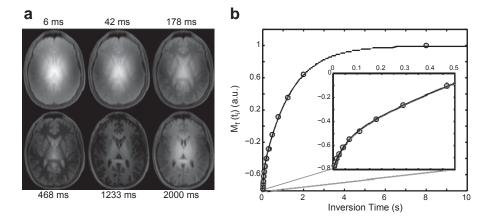


Fig. 5. Sample SIR-TFE images (a) and model fit (b) from a slice at the level of the lateral ventricles in a healthy control. (a) Images from six of the 14 inversion times are shown. Note the characteristic center brightening of the images due to Bi⁺ inhomogeneities. (b) Corresponding SIR data from a voxel in the genu of the corpus callosum. Note the agreement between the SIR data (circles) and biexponential model [solid black line, Eq. (7)] and the deviation from a monoexponential model, which is apparent at the shortest inversion times shown in the zoomed inset.

shorter t_d values were nearly identical to those at full recovery ($t_d = 10 \text{ s}$). Without this train, small deviations in PSR were observed at $t_d = 2.5 \text{ s}$; and these were more pronounced at $t_d = 1.25 \text{ s}$. Thus, the post-TFE saturation train allows for reduction of t_d (and scan times) with minimal parameter bias.

Representative 3D SIR-TFE data are shown in Fig. 5. Fig. 5a shows a single slice at the level of the lateral ventricles acquired at six of 14 t_i values. Note the characteristic center brightening due to B_{\perp}^+ inhomogeneities. Fig. 5b shows data from a single voxel in the genu of the corpus callosum and the corresponding model fit. Note the agreement between the SIR-TFE data and the biexponential model described by Eq. (7). Additionally, note the deviation from monoexponential recovery, which is especially evident at the shortest inversion times.

Based upon these fits, maps of qMT parameters were generated. Recall that these maps were filtered to reduce the impact of outliers. Fig. 6 displays representative qMT parameter maps without filtering, with the previously described locally-adaptive Gaussian filter, and with a global Gaussian filter. The locally-adaptive and global filters both removed outliers in the parameter maps (see arrow in the top row and the masks in the bottom row); however, the locally filtered maps were blurred to a much smaller degree. As a result, we employed the locally-adaptive approach herein. For all parameter maps, 14% of all voxels in the post-brain-extraction volume were smoothed using this approach. However, as seen in the bottom row of Fig. 6, a majority of these voxels were located along the brain surface or within the CSF.

Fig. 7 displays results from four representative slices in one healthy subject. The qMT parameters were uniform over most of the volume despite the presence of large $!1B_0\,\text{and/or}\,\,B_1^*$ field inhomogeneities (as indicated by the heterogeneity in the S_m maps). There does, however, appear to be some bias in the qMT parameter values in midbrain slices (black arrow), which typically (Moore et al., 2010) exhibit the largest field inhomogeneities and the lowest SNR. Nevertheless, these data suggest that robust qMT parameter mapping can be achieved throughout most of the brain using the 3D SIR-TFE protocol described herein.

ROIs were defined in a number of WM and GM regions as shown in Fig. 8. The boxplots in the top row of Fig. 9 display the mean ROI qMT parameters over the 12 healthy volunteers. For PSR, the mean value across all WM ROIs (17.6 \pm 1.3%) was higher than the values across all GM ROIs (10.3 \pm 1.6%). Additionally, heterogeneity within WM PSR values was observed, but should be interpreted with caution due to the effect of multiple comparisons. Nevertheless, differences between the following regions were detected: i) the genu of the corpus callosum and occipital WM (p = 0.026), ii) the genu of the corpus callosum and the corona radiata (p=0.026), iii) frontal and occipital

WM (p=0.041), and iv) frontal WM and the corona radiata (p=0.041). Fit k_{mf} values were higher in GM (24.4±4.4 s⁻¹) than in WM (14.5±1.5 s⁻¹). Additionally, note the large, biased k_{mf} values in and

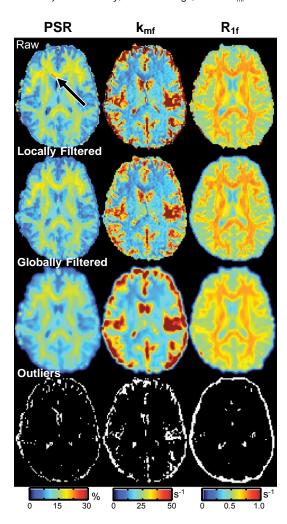


Fig. 6. Representative qMT parameter (PSR, k_{mf} , R_{1f}) maps with and without filtering. Shown are (1st row) raw parameter maps, (2nd row) parameter maps filtered with the locally-adaptive Gaussian filter, (3rd row) parameter maps filtered with the global Gaussian filter (with an identical kernel), and (4th row) masks of the outliers detected using the locally-adaptive filter. The arrow identifies a region with biased PSR values that are corrected by filtering. Note that the color-scale in these maps was chosen to highlight the outliers and is different than in Figs. 4 and 8.

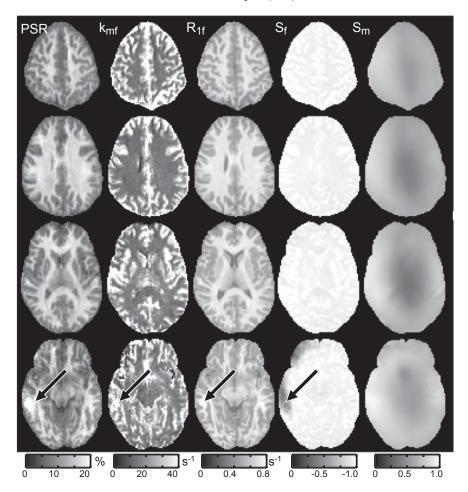


Fig. 7. Representative parameter maps from one subject (four of 30 slices are shown). The qMT parameters (PSR, k_{mf} , R_{1f}) and the inversion efficiency S_f were uniform over most of the volume despite the presence of large field inhomogeneities. There does, however, appear to be some bias in the qMT parameter values in midbrain slices (black arrows), which twoically exhibit the largest B_a and B_a^+ inhomogeneities. This results in a deviation of S_f from -1 in these regions.

around areas containing CSF, which is likely a consequence of the weak dependence of the signal on $k_{\rm mf}$ when PSR ≈ 0 (see the simulated data in Fig. 3c). For R_{1f} differences between WM (0.73 \pm 0.03 s $^{-1}$) and GM (0.58 \pm 0.05 s $^{-1}$) values were also observed. The boxplots in Fig. 9 give an indication of the variability of each parameter across the healthy cohort. To quantify this, the coefficient of variation was tabulated for each ROI, and the mean value across all ROI is given in Table 1. From this, it can be seen that the mean CV_cohort was b 10% for all of the qMT parameters, which is not surprising given the small age range of the healthy cohort scanned herein.

BA plots of the observed difference in mean ROI qMT parameters between scans are shown in the bottom row of Fig. 9; and the results from this analysis are given numerically in Table 1. The mean

difference for all ROIs across scans was close to zero for PSR (0.0%), $k_{\rm mf}$ (1.2 s $^{-1}$), and R_{1f} (0.01 s $^{-1}$), indicating a lack of bias and reasonable reproducibility. To further test this, a Wilcoxon signed-rank test was performed on the test–retest parameter values in each ROI. At the p = 0.05 level, no significant difference was observed between test and retest qMT parameters in any of the ROIs except for $k_{\rm mf}$ in the genu of the corpus callosum (p = 0.031). The test–retest coefficient of variation (CV $_{\rm retest}$) was also tabulated for each metric to further assess each parameter's variability across time. As shown in Table 1, the relative CV $_{\rm retest}$ values were consistent with the corresponding CV $_{\rm cohort}$ values, with $k_{\rm mf}$ exhibiting the highest variability. In terms of absolute CV values, the test–retest variability was approximately 20% lower than the across cohort variability.

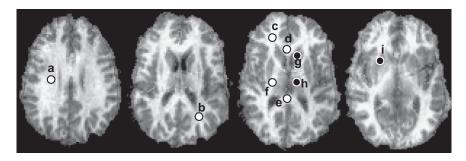


Fig. 8. Representative PSR maps from a single volunteer with corresponding ROIs (a = corona radiata; b = occipital WM; c = frontal WM; d = corpus callosum, genu; e = corpus callosum, splenium; f = internal capsule; g = head of caudate; h = thalamus; i = putamen). White and black dots represent WM and GM ROIs, respectively. In practice, ROIs were defined bilaterally and results were averaged across hemispheres. Here we show ROIs in one hemisphere for display purposes.

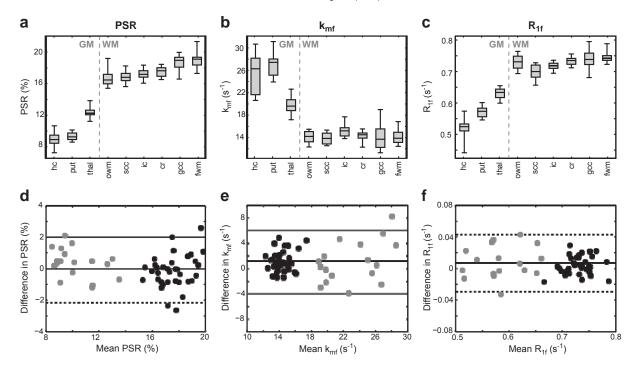


Fig. 9. (a-c) Boxplot of the mean ROI qMT parameters (hc = head of caudate; put = putamen; thal = thalamus; owm = occipital WM; scc = corpus callosum, splenium; ic = internal capsule; cr = corona radiata; gcc = corpus callosum, genu; fwm = frontal WM). On each box, the central mark is the median, the edges of the box are the 25th and 75th percentiles, and the whiskers extend to the most extreme data points. (d-f) Bland-Altman plots of the difference in parameters for WM (black) and GM (gray) ROIs across scans. The solid line is the mean difference, and the dashed lines are the limits of agreement (mean difference ±1.96 SD).

Discussion

This study demonstrates the feasibility of performing whole-brain gMT measurements in the human brain in vivo at high field. Pulsed saturation and SSFP-based approaches are difficult to implement at high field due to RF power limitations and/or magnetic field (B₁+ and !1Bo) inhomogeneities. In this study, we employed the SIR qMT approach, which has been suggested to be less sensitive to these issues. The biggest obstacles to overcome were i) the effect of B₁⁺ and !1Bo inhomogeneities on the inversion pulse and the readout and ii) the long scan times associated with SIR imaging. The former of these was mitigated by developing a novel B+ and !1B insensitive inversion composite pulse (Fig. 2) and employing a low-flip angle TFE readout; the latter was mitigated by the efficiency of the TFE readout along with additional protocol optimization (e.g., reducing the number of t_i values to 14, applying SENSE acceleration in two directions). Together this resulted in a robust (Fig. 9), whole-brain qMT imaging protocol with a scantime of less than 20 min.

Previous qMT imaging studies at lower field strengths (Dortch et al., 2011; Garcia et al., 2010; Gloor et al., 2008; Ropele et al., 2003; Sled and Pike, 2001; Sled et al., 2004; Yarnykh and Yuan, 2004) have reported PSR values in the range of 11–16% and 5–9% for WM and GM structures, respectively [PSR = F using the notation of Sled and Pike (2000, 2001) and M_{0b} using the notation of Henkelman et al. (1993)]. The PSR values presented herein (WM: 15–20%, GM: 9–

13%) were approximately 25% higher. PSR should be independent of field strength, so these differences may be related to the SIR-TFE sequence. As previously discussed, we modified the inversion pulse and readout of our 3.0-T SIR-FSE sequence to perform qMT imaging at 7.0 T. The effect of the TFE readout on the fit qMT parameters was assessed via numerical simulations and was found to result in little bias in PSR (Fig. 3). However, it should be noted that our previous report at 3.0 T employed a much longer TE (74 ms) than was employed herein (1.4 ms). Previous work (Bjarnason et al., 2005; Stanisz et al., 1999) has demonstrated that MT contrast is TE-dependent in WM due to the microanatomical compartition of water into myelin and nonmyelin water spaces. As a result, it is reasonable to assume that PSR may also exhibit a TE-dependence. In terms of the inversion pulse, we recognize that PSR is sensitive to the macromolecular pool lineshape and T_{2m} assumptions used in the numerical estimation of S_m. Similar to our previous studies (Dortch et al., 2011; Gochberg and Gore, 2007), we modeled the macromolecular pool using a Gaussian lineshape (T_{2m} = 10 ls) because the Super-Lorentzian exhibits an on-resonance singularity. Previous work using a 1-ms block inversion pulse at 3.0 T (Dortch et al., 2011) found that this was a reasonable approximation; however, this may not be true for the longer (5.5 ms), higher power composite inversion pulse employed herein. Additional work is needed to explore the field- and TE-dependence of PSR values obtained via the SIR technique. Nevertheless, the reported regional variation in PSR values was consistent with previous qMT imaging studies (Dortch et al., 2011;

Table 1 Test-retest reproducibility analysis of each qMT parameter (PSR, R_{1f} , and k_{mf}). Shown are the mean \pm SD parameter values across all ROIs for the test and retest scans, the resulting mean paired-difference between time-points, the limits-of-agreement (LOA), and the mean \pm SD test-retest coefficient of variation (CV_{retest}) across all ROIs. For comparison, the corresponding across-cohort coefficient of variation (CV_{cohort}) is also given.

| Parameter | Test scan (mean ± SD) | Retest scan (mean ± SD) | Difference | LOA | CV _{retest} (%) | CV _{cohort} (%) |
|------------------------------------|-----------------------|-------------------------|------------|----------------|--------------------------|--------------------------|
| PSR (%) | 15.2 ± 3.9 | 15.2 ±3.7 | 0.0 | (- 2.2, 2.1) | 4.9 ± 1.5 | 5.6 ± 1.9 |
| k _{mf} (s ⁻¹) | 16.8 ± 4.7 | 18.0 ±5.5 | 1.2 | (- 3.9, 6.3) | 8.2 ± 2.4 | 9.4 ± 3.6 |
| R _{1f} (s ⁻¹) | 0.68 ± 0.08 | 0.68 ±0.08 | 0.01 | (- 0.03, 0.04) | 1.9 ± 1.4 | 3.2 ± 1.3 |

Garcia et al., 2010; Sled et al., 2004; Underhill et al., 2009); and additional SIR-TFE studies in bovine serum albumin phantoms at 7.0 T (data not shown) found a linear relationship between macromolecular content and PSR. Thus, we postulate that the regional differences in PSR values reported herein are driven primarily by regional differences in myelin content, although the absolute values may be systematically larger than reported by other techniques.

Previous pulsed saturation and SSFP-based studies (Garcia et al., 2010; Gloor et al., 2008; Ropele et al., 2003; Sled and Pike, 2001; Sled et al., 2004; Yarnykh and Yuan, 2004) report k_{mf} values $[k_{mf} =$ $k_{\rm f}/\,F$ using the notation of Sled and Pike (2000, 2001); $k_{\rm mf}=R$ when $M_{of} = 1$ using the notation of Henkelman et al. (1993)] in the range of 20-40 s⁻¹ across the brain. A previous SIR study (Dortch et al., 2011) at 3.0 T reports k_{mf} values that are approximately 2-fold slower (10-15 s⁻¹) with values that are slower in WM than GM, which is consistent with the results presented herein. The discrepancies between techniques are not surprising given the reported difficulty of using pulsed saturation to determine k_{mf} (Portnoy and Stanisz, 2007). In terms of the current study, it should be noted that k_{mf} showed the largest variability of the qMT parameters, which is consistent with the results from the Monte Carlo simulations. We do not expect this to be a significant drawback as k_{mf} has been shown to be insensitive to the pathological changes in spinal cord WM (Smith et al., 2009).

While there have been no previous reports of R_{1f} in human brain at 7.0 T, it can be shown that the observed T_1 typically reported is $\approx 1/R_{1f}$. Using this relationship, the mean WM and GM observed T_1 values were 1372 and 1724 ms, respectively, which are within the range of previously reported values in human brain at 7.0 T (Wright et al., 2008). As expected, we noted a significant correlation between R_{1f} and PSR in the healthy human brain (data not shown); however, T_1 is sensitive to overall tissue composition [e.g., water content (Kiricuta and Simplaceanu, 1975)] and is believed to be a less specific marker for myelin in WM.

The increased SNR available at 7.0 T was used here to decrease scan time (\approx 40 s/slice) and increase resolution (2 × 2 × 3 mm³) relative to our 3.0-T protocol. Moving forward, it may be advantageous to look at higher resolution protocols. If we assume that all imaging parameters are the same, increasing the resolution to 1×1×3 mm³ would result in an approximately two-fold decrease in SNR (\approx 70 at thermal equilibrium, assuming we increase the number of acquired points to hold the field-of-view constant). Based upon previous simulation work (Li et al., 2010) as well as the simulation work presented herein, this would be sufficient to robustly fit qMT parameters over most of the brain. Thus, it appears that high-resolution qMT imaging may be feasible in the human brain in vivo at 7.0 T using a protocol similar to that described herein.

Conclusions

The results of this study demonstrate the feasibility of performing qMT imaging in human brain in vivo at high field. The developed SIR-TFE protocol allowed for whole-brain qMT imaging in less than 20 min. In healthy subjects, intra-subject reliability (i.e., test–retest) was demonstrated despite large $!1B_0\, {\rm and}\,\, B_l^+$ variations. Additionally, a high level of inter-subject reproducibility was demonstrated for the qMT parameters. Future work includes investigating high-resolution protocols to look at cortical features of qMT parameters and application of the approach in a cohort of multiple sclerosis patients.

Acknowledgments

This work was supported by NIH/NBIB K01 EB009120 (SAS), NIH T32 EB001628 (JCG), NIH EB00461 (JCG), and Vanderbilt Bridge Funding (DFG).

References

- Berry, I., Barker, G., Barkhof, F., Campi, A., Dousset, V., Franconi, J., Gass, A., Schreiber, W., Miller, D., Tofts, P., 1999. A multicenter measurement of magnetization transfer ratio in normal white matter. J. Magn. Reson. Imaging 9, 441–446.
- Bjarnason, T., Vavasour, I., Chia, C., Mackay, A., 2005. Characterization of the NMR behavior of white matter in bovine brain. Magn. Reson. Med. 54, 1072–1081.
- Bruno, S.D., Barker, G.J., Cercignani, M., Symms, M., Ron, M.A., 2004. A study of bipolar disorder using magnetization transfer imaging and voxel-based morphometry. Brain 127, 2433–2440.
- Catalaa, I., Grossman, R.I., Kolson, D.L., Udupa, J.K., Nyul, L.G., Wei, L., Zhang, X., Polansky, M., Mannon, L.J., McGowan, J.C., 2000. Multiple sclerosis: magnetization transfer histogram analysis of segmented normal-appearing white matter. Radiology 216, 351–355.
- Cercignani, M., Alexander, D., 2006. Optimal acquisition schemes for in vivo quantitative magnetization transfer MRI. Magn. Reson. Med. 56, 803–810.
- Constable, R.T., Gore, J.C., 1992. The loss of small objects in variable TE imaging: implications for FSE, RARE, and EPI. Magn. Reson. Med. 28, 9–24.
- Dortch, R.D., Li, K., Gochberg, D.F., Welch, E.B., Dula, A.N., Tamhane, A.A., Gore, J.C. Smith, S.A., 2011. Quantitative magnetization transfer imaging in human brain at 3 T via selective inversion recovery. Magn. Reson. Med. 66, 1346–1352.
- Dousset, V., Grossman, R.I., Ramer, K.N., Schnall, M.D., Young, L.H., Gonzalez-Scarano, F., Lavi, E., Cohen, J.A., 1992. Experimental allergic encephalomyelitis and multiple sclerosis: lesion characterization with magnetization transfer imaging. Radiology 182. 483–491.
- Edzes, H.T., Samulski, E.T., 1977. Cross relaxation and spin diffusion in the proton NMR or hydrated collagen. Nature 265, 521–523.
- Filippi, M., Rocca, M.A., 2004. Magnetization transfer magnetic resonance imaging in the assessment of neurological diseases. J. Neuroimaging 14, 303–313.
- Garcia, M., Gloor, M., Wetzel, S.G., Radue, E.-W., Scheffler, K., Bieri, O., 2010. Characterization of normal appearing brain structures using high-resolution quantitative magnetization transfer steady-state free precession imaging. NeuroImage 52, 532–537.
- Gass, A., Barker, G.J., Kidd, D., Thorpe, J.W., MacManus, D., Brennan, A., Tofts, P.S., Thompson, A.J., McDonald, W.I., Miller, D.H., 1994. Correlation of magnetization transfer ratio with clinical disability in multiple sclerosis. Ann. Neurol. 36, 62–67.
- Gloor, M., Scheffler, K., Bieri, O., 2008. Quantitative magnetization transfer imaging using balanced SSFP. Magn. Reson. Med. 60, 691–700.
- Gochberg, D., Gore, J., 2007. Quantitative magnetization transfer imaging via selective inversion recovery with short repetition times. Magn. Reson. Med. 57, 437–441.
- Gochberg, D., Kennan, R., Gore, J., 1997. Quantitative studies of magnetization transfer by selective excitation and T_1 recovery. Magn. Reson. Med. 38, 224–231.
- Graham, S.J., Henkelman, R.M., 1997. Understanding pulsed magnetization transfer. J. Magn. Reson. Imaging 7, 903–912.
- Henkelman, R., Huang, X., Xiang, Q., Stanisz, G., Swanson, S., Bronskill, M., 1993. Quantitative interpretation of magnetization transfer. Magn. Reson. Med. 29, 759–766.
- Kabani, N.J., Sled, J.G., Chertkow, H., 2002a. Magnetization transfer ratio in mild cognitive impairment and dementia of Alzheimer's type. NeuroImage 15, 604–610. Kabani, N.J., Sled, J.G., Shuper, A., Chertkow, H., 2002b. Regional magnetization transfer ratio change in mild cognitive impairment. Magn. Recognitive impairment.
 - transfer ratio changes in mild cognitive impairment. Magn. Reson. Med. 47, 143–148.
- Kalkers, N.F., Hintzen, R.Q., van Waesberghe, J.H., Lazeron, R.H., van Schijndel, R.A., Ader, H.J., Polman, C.H., Barkhof, F., 2001. Magnetization transfer histogram parameters reflect all dimensions of MS pathology, including atrophy. J. Neurol. Sci. 184, 155–162.
- Kiricuta Jr., I.C., Simplaceanu, V., 1975. Tissue water content and nuclear magnetic resonance in normal and tumor tissues. Cancer Res. 35, 1164–1167.
- Koenig, S.H., 1991. Cholesterol of myelin is the determinant of gray—white contrast in MRI of brain. Magn. Reson. Med. 20, 285–291.
- Kucharczyk, W., Macdonal, P., Stanisz, G., Henkelman, R., 1994. Relaxivity and magnetization transfer of white matter lipids at MR imaging: importance of cerebrosides and pH. Radiology 192, 521–529.
- Levesque, I.R., Sled, J.G., Pike, G.B., 2011. Iterative optimization method for design of quantitative magnetization transfer imaging experiments. Magn. Reson. Med. 66, 635–643.
- Li, K., Zu, Z., Xu, J., Janve, V.A., Gore, J.C., Does, M.D., Gochberg, D.F., 2010. Optimized inversion recovery sequences for quantitative T₁ and magnetization transfer imaging. Magn. Reson. Med. 64, 491–500.
- Moore, J., Jankiewicz, M., Zeng, H., Anderson, A.W., Gore, J.C., 2010. Composite RF pulses for B_1^+ -insensitive volume excitation at 7 Tesla. J. Magn. Reson. 205, 50–62.
- Morrison, C., Stanisz, G., Henkelman, R.M., 1995. Modeling magnetization transfer for biological-like systems using a semi-solid pool with a super-Lorentzian lineshape and dipolar reservoir. J. Magn. Reson. B 108, 103–113.
- Mougin, O.E., Coxon, R.C., Pitiot, A., Gowland, P.A., 2010. Magnetization transfer phenomenon in the human brain at 7 T. Neurolmage 49, 272–281.
- Mugler III, J.P., Brookeman, J.R., 1990. Three-dimensional magnetization-prepared rapid gradient-echo imaging (3D MP RAGE). Magn. Reson. Med. 15, 152–157.
- Odrobina, E., Lam, T., Pun, T., Midha, R., Stanisz, G., 2005. MR properties of excised neural tissue following experimentally induced clemyelination. NMR Biomed. 18, 277–284.
- Ou, X., Sun, S.W., Liang, H.F., Song, S.K., Gochberg, D.F., 2009. The MT pool size ratio and the DTI radial diffusivity may reflect the myelination in shiverer and control mice. NMR Biomed. 22, 480–487.
- Pike, G.B., 1996. Pulsed magnetization transfer contrast in gradient echo imaging: a two-pool analytic description of signal response. Magn. Reson. Med. 36, 95–103.
- Portnoy, S., Stanisz, G.J., 2007. Modeling pulsed magnetization transfer. Magn. Reson. Med. 58, 144–155.

- Ramani, A., Dalton, C., Miller, D.H., Tofts, P.S., Barker, G.J., 2002. Precise estimate of fundamental in-vivo MT parameters in human brain in clinically feasible times. Magn. Reson. Imaging 20, 721–731.
- Ropele, S., Seifert, T., Enzinger, C., Fazekas, F., 2003. Method for quantitative imaging of the macromolecular ¹H fraction in tissues. Magn. Reson. Med. 49, 864–871.
- Schmierer, K., Scaravilli, F., Altmann, D., Barker, G., Miller, D., 2004. Magnetization transfer ratio and myelin in postmortem multiple sclerosis brain. Ann. Neurol. 56, 407–415.
- Schmierer, K., Tozer, D.J., Scaravilli, F., Altmann, D.R., Barker, G.J., Tofts, P.S., Miller, D.H., 2007. Quantitative magnetization transfer imaging in postmortem multiple sclerosis brain. J. Magn. Reson. Imaging 26, 41–51.
- Sled, J.G., Pike, G.B., 2000. Quantitative interpretation of magnetization transfer in spoiled gradient echo MRI sequences. J. Magn. Reson. 145, 24–36.
- Sled, J.G., Pike, G.B., 2001. Quantitative imaging of magnetization transfer exchange and relaxation properties in vivo using MRI. Magn. Reson. Med. 46, 923–931.
- Sled, J.G., Levesque, I., Santos, A.C., Francis, S.J., Narayanan, S., Brass, S.D., Arnold, D.L., Pike, G.B., 2004. Regional variations in normal brain shown by quantitative magnetization transfer imaging. Magn. Reson. Med. 51, 299–303.
- Smith, S.M., 2002. Fast robust automated brain extraction. Hum. Brain Mapp. 17, 143–155.
- Smith, S., Golay, X., Fatemi, A., Mahmood, A., Raymond, G., Moser, H., van Zijl, P., Stanisz, G., 2009. Quantitative magnetization transfer characteristics of the

- human cervical spinal cord in vivo: application to adrenomyeloneuropathy. Magn. Reson. Med. 61, 22–27.
- Stanisz, G., Kecojevic, A., Bronskill, M., Henkelman, R., 1999. Characterizing white matter with magnetization transfer and T_2 . Magn. Reson. Med. 42, 1128–1136.
- Underhill, H.R., Yuan, C., Yarnykh, V.L., 2009. Direct quantitative comparison between cross-relaxation imaging and diffusion tensor imaging of the human brain at 3.0 T. Neurolmage 47, 1568–1578.
- Underhill, H.R., Rostomily, R.C., Mikheev, A.M., Yuan, C., Yarnykh, V.L., 2011. Fast bound pool fraction imaging of the in vivo rat brain: association with myelin content and validation in the C6 glioma model. NeuroImage 54, 2052–2065.
- Viola, P., Wells, W.M., 1997. Alignment by maximization of mutual information. Int. J. Comput. Vis. 24, 137–154.
- Wolff, S.D., Balaban, R.S., 1989. Magnetization transfer contrast (MTC) and tissue water proton relaxation in vivo. Magn. Reson. Med. 10, 135–144.
- Wright, P.J., Mougin, O.E., Totman, J.J., Peters, A.M., Brookes, M.J., Coxon, R., Morris, P.E., Clemence, M., Francis, S.T., Bowtell, R.W., Gowland, P.A., 2008. Water proton T₁ measurements in brain tissue at 7, 3, and 1.5 T using IR-EPI, IR-TSE, and MPRAGE: results and optimization. MAGMA 21, 121–130.
- Yarnykh, V.L., 2007. Actual flip-angle imaging in the pulsed steady state: a method for rapid three-dimensional mapping of the transmitted radiofrequency field. Magn. Reson. Med. 57, 192–200.
- Yarnykh, V.L., Yuan, C., 2004. Cross-relaxation imaging reveals detailed anatomy of white matter fiber tracts in the human brain. NeuroImage 23, 409–424.